



**Thermo Scientific**

# **Acclaim Carbonyl C18**

## **Column Product Manual**

**P/N: 065451-02      April 2013**

## **Product Manual**

**for**

### **Acclaim Carbonyl C18 RSLC, Analytical Columns**

2.2  $\mu\text{m}$ , 2.1 x 100mm (P/N 077972)

2.2  $\mu\text{m}$ , 2.1 x 150mm (P/N 077973)

2.2  $\mu\text{m}$ , 3.0 x 100mm (P/N 077974)

### **Acclaim Carbonyl C18, Analytical Columns**

5  $\mu\text{m}$ , 4.6 x 250 mm (P/N083214)

5  $\mu\text{m}$ , 4.6 x 150 mm (P/N 079008)

3  $\mu\text{m}$ , 3.0 x 250mm (P/N 079009)

3  $\mu\text{m}$ , 3.0 x 150mm (P/N 079010)

3  $\mu\text{m}$ , 2.1 x 150mm (P/N 079011)

### **Acclaim Carbonyl C18, Guard Columns**

5  $\mu\text{m}$ , 4.6 x 10mm (P/N 079014)

5  $\mu\text{m}$ , 3.0 x 10mm (P/N 079013)

5  $\mu\text{m}$ , 2.1 x 10mm (P/N 079012)

© 2013 Thermo Fisher Scientific Inc. All rights reserved.

All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries.

Thermo Fisher Scientific Inc. provides this document to its customers with a product purchase to use in the product operation. This document is copyright protected and any reproduction of the whole or any part of this document is strictly prohibited, except with the written authorization of Thermo Fisher Scientific Inc.

The contents of this document are subject to change without notice. All technical information in this document is for reference purposes only. System configurations and specifications in this document supersede all previous information received by the purchaser.

Thermo Fisher Scientific Inc. makes no representations that this document is complete, accurate or error free and assumes no responsibility and will not be liable for any errors, omissions, damage or loss that might result from any use of this document, even if the information in the document is followed properly.

This document is not part of any sales contract between Thermo Fisher Scientific Inc. and a purchaser. This document shall in no way govern or modify any Terms and Conditions of Sale, which Terms and Conditions of Sale shall govern all conflicting information between the two documents.

***This document is for Research Use Only. This document is Not for use in Diagnostic Procedures.***

Revision History:

Revision 02, April, 2013, Reformatted for Thermo Scientific. Added 3 $\mu$ m and 5 $\mu$ m columns to the product line.

## Safety and Special Notices

Make sure you follow the precautionary statements presented in this guide. The safety and other special notices appear in boxes.

Safety and special notices include the following:



**SAFETY**

*Indicates a potentially hazardous situation which, if not avoided, could result in death or serious injury.*



**WARNING**

*Indicates a potentially hazardous situation which, if not avoided, could result in damage to equipment.*



**CAUTION**

*Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury. Also used to identify a situation or practice that may seriously damage the instrument, but will not cause injury.*



**NOTE**

*Indicates information of general interest.*

**IMPORTANT**

*Highlights information necessary to prevent damage to software, loss of data, or invalid test results; or might contain information that is critical for optimal performance of the system.*

**Tip**

*Highlights helpful information that can make a task easier.*

# Contents

<b>1. Introduction</b>	<b>6</b>
<b>2. Installation</b>	<b>7</b>
2.1 Preparation of the Mobile Phase	7
2.1.1 Ammonium Acetate buffer	7
2.2 HPLC System Set-up	7
2.3 Conditioning the Column	7
2.4 Ensuring Column Performance	7
2.5 Optimizing Chromatographic Conditions	8
2.5.1 Standard Acetonitrile conditions	9
2.5.2 Standard Methanol conditions	10
2.6 Real Sample Analysis	11
2.7 Column Storage	11
<b>3. Column Care</b>	<b>12</b>
3.1 General Guidelines	12
3.1.1 Recommended Ranges of Operation	12
3.2 Recommended Operating pH Range	12
3.3 Recommended Operating Temperature	12
3.4 Recommended Flow Rate	12
3.5 Column Washing Procedure	14
<b>4. Frequently Asked Questions (FAQ)</b>	<b>15</b>
<b>5. Applications</b>	<b>17</b>

# 1. Introduction

Aldehydes and ketones are common pollutants in air and water. The analytical difficulties that need to be overcome include their volatility, their reactivity, and their modest UV absorption. The reaction with dinitrophenylhydrazine (DNPH) is a convenient means of trapping, stabilizing, and tagging these substances. Several standard methods have been developed to apply this chemistry to various environmental situations. Some of the better known ones include CARB 1004 for vehicle exhaust, EPA 554 for drinking water, EPA 1667 for pharmaceutical wastewater, and EPA 8315 for general wastewater. Each standard method has a particular target compound list.

The Thermo Scientific™ Acclaim™ Carbonyl columns are silica-based reversed phase columns designed specifically for separating DNPH derivatives of aldehydes and ketones. They exhibit superior resolution compared with other commercially available columns. All the popular standard methods can be performed with the Acclaim Carbonyl columns.

The Acclaim Carbonyl columns are available in three particle sizes, 5 µm, 3 µm, and 2.2 µm. The 5 µm, 4.6 x 250 mm column is intended for traditional HPLC equipment, and provides excellent resolution. The 3 µm, 3.0x150 mm column gives faster run times and lower solvent consumption, and is compatible with most traditional HPLC systems. The 2.2 µm columns are designed for newer UHPLC systems and feature excellent resolution, speed and solvent consumption.

Acclaim Carbonyl C18 Analytical Columns		
Particle Size	Dimensions	Part Number
5 µm	4.6 x 250 mm	083214
	4.6 x 150 mm	079008
	3.0 x 250 mm	079009
3 µm	3.0 x 150 mm	079010
	2.1 x 150 mm	079011
	2.1 x 100 mm	077972
2.2 µm	2.1 x 150 mm	077973
	3.0 x 100 mm	077974

Acclaim Carbonyl C18 Guard Columns	
Dimensions	Part Number
4.6 x 10 mm	079014
3.0 x 10 mm	079013
2.1 x 10 mm	079012
Holder, V2	069580

## 2. Installation

### 2.1 Preparation of the Mobile Phase

The recommended mobile phase system contains an aqueous portion and a solvent portion. The aqueous portion can be D.I. water or a buffer (e.g., 2 mM ammonium acetate buffer). While D.I. water works well in most cases, a buffered aqueous solution sometimes may yield slightly better consistency of results. The solvent can be acetonitrile or methanol, depending on the application. All ingredients for the mobile phases are available in “HPLC grade” or better and should be used to make mobile phases. If you use an in-house water purifier, be sure it is correctly maintained.

#### 2.1.1 Preparation of 2 mM ammonium acetate buffer

1. Measure  $57 \pm 3 \mu\text{L}$  of glacial acetic acid and  $154 \pm 5 \text{ mg}$  of ammonium acetate into  $1000 \pm 5 \text{ mL}$  of deionized water.

### 2.2 HPLC System Set-up

The Acclaim Carbonyl columns may be used on any standard HPLC system equipped with an HPLC pump, a column oven, a UV detector, and an injector. To minimize ‘downtime’, ensure that the whole system is primed before starting your column conditioning.

### 2.3 Conditioning the Column

When installing a new Carbonyl column for the first time, wash the column with pure methanol for ~ 20 column volumes (or 50 minutes at 1 mL/min), and send the effluent directly to waste. Reconnect the column to the detector, and equilibrate it with the desired mobile phase for at least 20 column volumes before making your first injection.



NOTE

*Thermo Scientific recommends that you always read the manual for a new column before installing it for the first time. The manual contains information regarding the operational limits of the column, as well as advice on how to optimize your separation.*

### 2.4 Ensuring Column Performance

Before running any samples, Thermo Scientific recommends that you first confirm the performance of the column by reproducing the lot validation report chromatogram shipped with the column. Compare your results with the one reported in the quality assurance report. At least three injections should be made.



NOTE

*Due to various reasons, such as differences in HPLC systems, mobile phases, oven temperature control, etc., you may observe somewhat different peak resolution from those shown in the report. In this case, please contact us with your test chromatogram for technical support and/or optimize chromatographic conditions using the methods stated in “Optimizing Chromatographic Conditions”.*

## 2.5 Optimizing Chromatographic Conditions

Some typical chromatograms using the following conditions are shown in sections 2.5.1 and 2.5.2.



**NOTE**

*Mobile phase composition and oven temperature are the two main factors influencing the separation.*



**NOTE**

*The mobile phase gradients listed below are intended to provide a starting point, and should be modified as necessary for optimal separation.*



## 2 – Installation

### 2.5.1 Standard Acetonitrile conditions

These gradients are satisfactory for CARB 1004, EPA 554, EPA 1667 and EPA 8315. Depending on the gradient delay volume of your system, you may want to adjust the times for the start and end of the gradient. You may also want to adjust the equilibration time between injections. For EPA 1667, the gradient is optional; you may use simply the initial isocratic conditions if there are no late-eluting interferences.

Column (Notes)	Flow Rate (mL/min)	Temp. (°C)	Data Rate (Hz)	Inj. Volume (µL)	Gradient		
					Time (min)	% Acetonitrile	% Aqueous
5 µm, 4.6 x 250 mm  Standard HPLC (Fig. 7)	1.50	28	5	10	-9	52	48
					0	52	48
					20	52	48
					28	90	10
					32	90	10
3 µm, 3.0 x 150 mm  Standard HPLC; optimized for speed (Fig. 5)	0.90	28	10	5	-3.5	53	47
					0.0	53	47
					7.5	53	47
					11.5	90	10
					16.0	90	10
3 µm, 3.0 x 250 mm  Standard HPLC; optimized for resolution	0.60	28	5	5	-12.0	53	47
					0.0	53	47
					20.5	53	47
					29.0	90	10
					34.0	90	10
3 µm, 2.1 x 150 mm  Micro HPLC; optimized for low solvent consumption	0.30	28	10	2	-6.0	53	47
					0.0	53	47
					11.0	53	47
					19.0	90	10
					24.0	90	10
2.2 µm, 2.1 x 150 mm Low gradient delay UHPLC; optimized for resolution (Fig. 4)	0.40	28	25	1	-4.5	52	48
					0.0	52	48
					8.3	52	48
					15.0	100	0
					18.0	100	0
2.2 µm, 2.1 x 100 mm Low gradient delay UHPLC; optimized for speed (Fig. 1 & 3)	0.75	28	25	1	-1.7	52	48
					0.0	52	48
					2.9	52	48
					5.3	100	0
					6.2	100	0
2.2 µm, 3.0 x 100 mm  Quaternary HPLC	1.00	28	25	2	-3.0	52	48
					0.0	52	48
					4.4	52	48
					8.0	100	0
					9.8	100	0

## 2.5.2 Standard Methanol conditions

These gradients are satisfactory for EPA 554. Depending on the gradient delay volume of your system, you may want to adjust the times for the start and end of the gradient. You may also want to adjust the equilibration time between injections.

Column (Notes)	Flow Rate (mL/min)	Temp. (°C)	Data Rate (Hz)	Inj. Volume ( $\mu$ L)	Gradient		
					Time (min)	% Methanol	% Aqueous
5 $\mu$ m, 4.6 x 250 mm  Standard HPLC	1.25	40	5	20	-9	72	28
					0	72	28
					20	72	28
					33	100	0
					40	100	0
3 $\mu$ m, 3.0 x 150 mm  Standard HPLC; optimized for speed (Fig. 6)	0.80	40	10	5	-6.0	72	28
					0.0	72	28
					7.0	72	28
					12.5	100	0
					16.0	100	0
3 $\mu$ m, 3.0 x 250 mm  Standard HPLC; high resolution	0.53	40	2	5	-12.0	72	28
					0.0	72	28
					19.0	72	28
					32.0	100	0
					40.0	100	0
3 $\mu$ m, 2.1 x 150 mm  Micro HPLC; low solvent consumption	0.26	40	5	1	-12.0	72	28
					0.0	72	28
					11.5	72	28
					14.5	100	0
					24.0	100	0
2.2 $\mu$ m, 2.1 x 150 mm  Low gradient delay UHPLC	0.50	42	25	1	-1.9	70	30
					0.0	70	30
					5.8	70	30
					8.9	100	0
					10.5	100	0
2.2 $\mu$ m, 2.1 x 100 mm  Low gradient delay UHPLC (Fig. 2)	0.50	42	25	1	-1.7	70	30
					0.0	70	30
					3.4	70	30
					5.5	100	0
					7.0	100	0
2.2 $\mu$ m, 3.0 x 100 mm  Quaternary HPLC	0.75	42	10	2	-3.5	70	30
					0.0	70	30
					3.5	70	30
					7.0	100	0
					9.3	100	0

## 2.6 Real Sample Analysis

Once satisfactory results have been obtained using your test mix, you are ready to run samples. The same conditions that separate the test mix should be used to analyze your samples.

## 2.7 Column Storage

After use, the column can be stored in the mobile phase for short periods of time (e.g. overnight). For longer term storage (longer than one week), it is recommended that you store the column in pure methanol or acetonitrile.

## 3. Column Care

### 3.1 General Guidelines

These columns should be used with the same precautions you would take for any other silica-based reversed phase columns. Please refer to the table below for recommended operational guidelines.

#### 3.1.1 Recommended Ranges of Operation

Particle Size	Column Dimensions	Maximum Pressure (bar)	Maximum Flow Rate (mL/min)	Typical Flow Rate (mL/min)	pH Range	Typical Temperature (° C)	Maximum Temperature (° C)
5 µm	4.6 x 250 mm	400	2.0	0.8 – 1.5	2.5 – 7.5	25 – 35	50
	4.6 x 150 mm	400	2.0	0.8 – 1.5	2.5 – 7.5	25 – 35	50
3 µm	3.0 x 250 mm	800	1.0	0.4 – 0.8	2.5 – 7.5	25 – 35	50
	3.0 x 150 mm	600	1.0	0.4 – 0.8	2.5 – 7.5	25 – 35	50
	2.1 x 150 mm	600	0.5	0.2 – 0.4	2.5 – 7.5	25 – 35	50
2.2 µm	2.1 x 150 mm	800	0.1	0.25 – 0.75	2.5 – 7.5	25 – 35	50
	2.1 x 100 mm	600	0.1	0.25 – 0.75	2.5 – 7.5	25 – 35	50
	3.0 x 100 mm	550	1.5	0.4 – 1.2	2.5 – 7.5	25 – 35	50

### 3.2 Recommended Operating pH Range

The pH of the mobile phase has little effect on retention times or selectivity on the Acclaim Carbonyl columns, and therefore need not be varied for most samples. To ensure the longest possible lifetime for these columns, a mobile phase that is 'silica friendly' should be used (pH 2.5 – 7.5). In most cases, a simple methanol (or acetonitrile)/water (or ammonium acetate) mobile phase system will work very well.

### 3.3 Recommended Operating Temperature

The separation of carbonyl DNPH is moderately sensitive to changes in temperature. We have found that optimal separation occurs between 25 °C and 35 °C for most applications using acetonitrile gradients, and on most systems. When using methanol gradients, moderately elevated temperatures (35 – 45 °C) are recommended to reduce the back pressure; this permits conventional HPLC equipment to use the 2.2 µm RSLC columns without exceeding the pressure limits of the system. Although the Acclaim Carbonyl columns can be used within a broader temperature range, we have found no practical reason to use them outside the recommended range in order to improve the separation. On the other hand, a mild operating temperature helps to prolong column lifetime.

### 3.4 Recommended Flow Rate

It is extremely important not to expose the columns to surges in column pressure. When starting up a system from idle, for a 4.6-mm i.d. column, gradually increase the flow rate from 0.5 mL/min up to the desired flow rate in 0.1-0.2 mL/min increments.



### 3.5 Column Washing Procedure

All samples should be pre-treated and filtered before being injected onto the column. In addition, a guard column is recommended for real sample analysis to prolong the lifetime of the analytical column. If the column does need to be cleaned, such as after long-term storage, the following procedure can be used as a guideline. The flow rates given here are for a 4.6 mm i.d. column; for a 3 mm i.d. column reduce the flow to 40% and for a 2.1 mm i.d. column reduce the flow to 20%.

1. Equilibrate the column with methanol/water v/v 50/50 for 10 column volumes at 0.5 to 1 mL/min.
2. Then, wash the column with pure methanol, acetonitrile, or acetone for 20 column volumes at 0.5 to 1 mL/min.
3. Finally, wash the column with methanol/water v/v 80/20 for 10 column volumes.
4. Before any injection is made, the column should be equilibrated with your initial mobile phase composition for at least 20 column volumes.

## 4. Frequently Asked Questions (FAQ)

1. How do the Acclaim Carbonyl columns compare with other columns for Aldehydes and Ketones?

The Acclaim Carbonyl column offers unique selectivity for DNPH derivatives. There is no need for excessively long columns or run times, or for complex gradients to compensate for inadequate selectivity. The Acclaim Carbonyl columns offer a combination of selectivity, resolution and speed that is unmatched by any other column.

2. Which particle size / format should I use?

The 5- $\mu$ m particle size is used for standard applications when methods are being developed or if the sample is complex. The 3- $\mu$ m particle size is used for faster, higher efficiency separations. The 2.2- $\mu$ m particle size is used in UHPLC systems that operate at higher pressure, and provide high-speed analysis. The 250-mm long columns are optimum for highest resolution and the 150-mm long columns are for high speed separations. The 2.1-mm internal diameter columns are the preferred column for LC/MS applications, or in cases where sample size is limited and greater sensitivity is required, or when solvent savings are desired.

3. My chromatogram shows bad peak shapes and low efficiencies. What is the problem, and how do I resolve it?

- A. Check the system, connections, and tubing for excessive extra column volume. Fix and replace if needed.
- B. Be sure the column is fully equilibrated with mobile phase.
- C. Make sure you are not injecting too large a sample or a sample in the wrong solvent.
- D. Run the column performance test described in the Quality Assurance Report (QAR). Replace the column if it is necessary.

4. Why am I observing high column backpressure?

- A. Check the injection valve and tubing for possible clogging.
- B. Wash the column using the protocol described in “Section 3.5.”
- C. Use a guard column and/or pre-column filter and replace it on a regular basis.

5. Why is the selectivity on my column different from the Quality Assurance Report (QAR), when using the condition described in QAR?

- A. Check your mobile phase composition.
- B. If you are proportioning, try using different lines, or pre-mix the solution and check the selectivity.
- C. Check your column temperature (oven temperature) and make sure it has been calibrated recently.
- D. Use the column within its pH limits. Failing to do so will result in undesirable selectivity change.

6. Should I use D.I. water or buffer for the analysis?

Both will give good result.. There is virtually no difference in results between D.I. water and buffer. We have observed that in some cases, a slightly buffered aqueous solution may yield marginally better consistency


7. When should I use methanol instead of acetonitrile?

The chemistry of DNPH derivatives produces isomers for some aldehydes and ketones. These isomers are better resolved with acetonitrile; with methanol, the isomers tend to coelute. (This is a unique feature of the Acclaim Carbonyl columns.) Some pairs of analytes are better resolved using acetonitrile, especially in CARB 1004 and EPA 8315. Acetonitrile usually permits moderately faster analysis times. Methanol is less expensive and less subject to supply problems than acetonitrile. Methanol is more viscous than acetonitrile, and for some column formats, may result in excess pressure. The reference method in EPA 554 uses methanol, and this may be the more comfortable option for regulated laboratories.




## 5. Applications


Figure 1 Rapid DNPH Aldehyde and Ketone Standards using Acclaim Carbonyl RSLC with Acetonitrile Mobile Phase

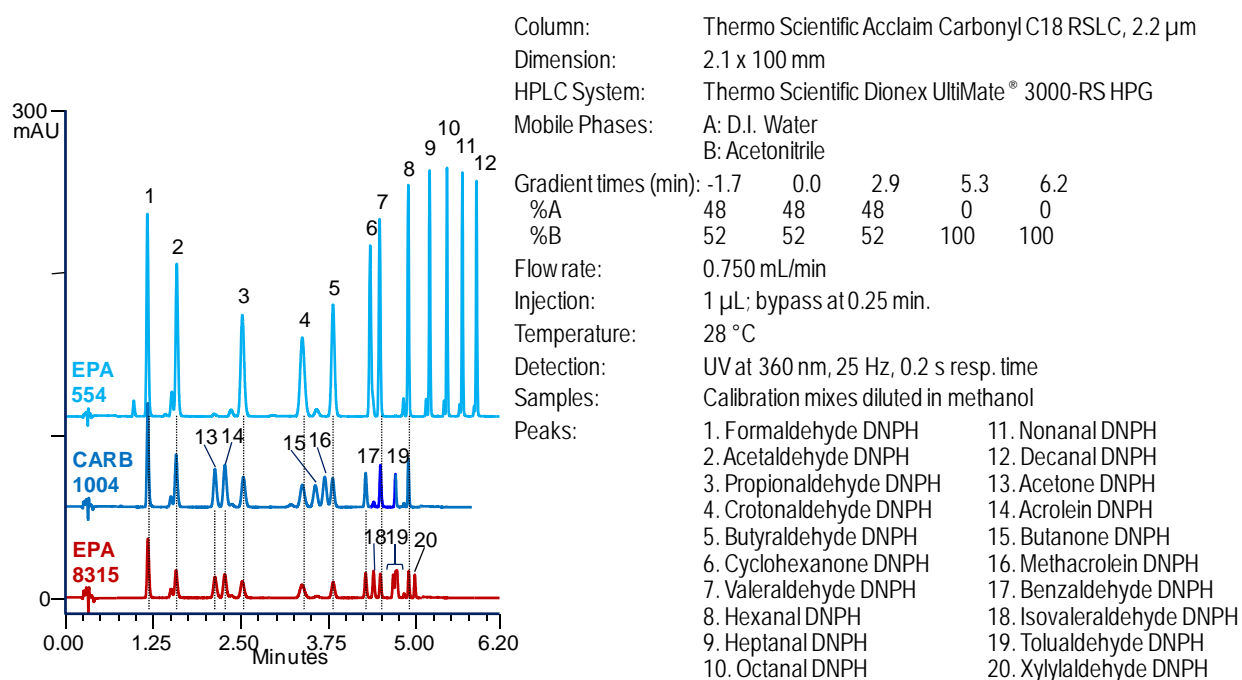


**Figure 2** Rapid EPA 554 DNPH Aldehyde and Ketone Standards using Acclaim Carbonyl RSLC with Acetonitrile Mobile Phase

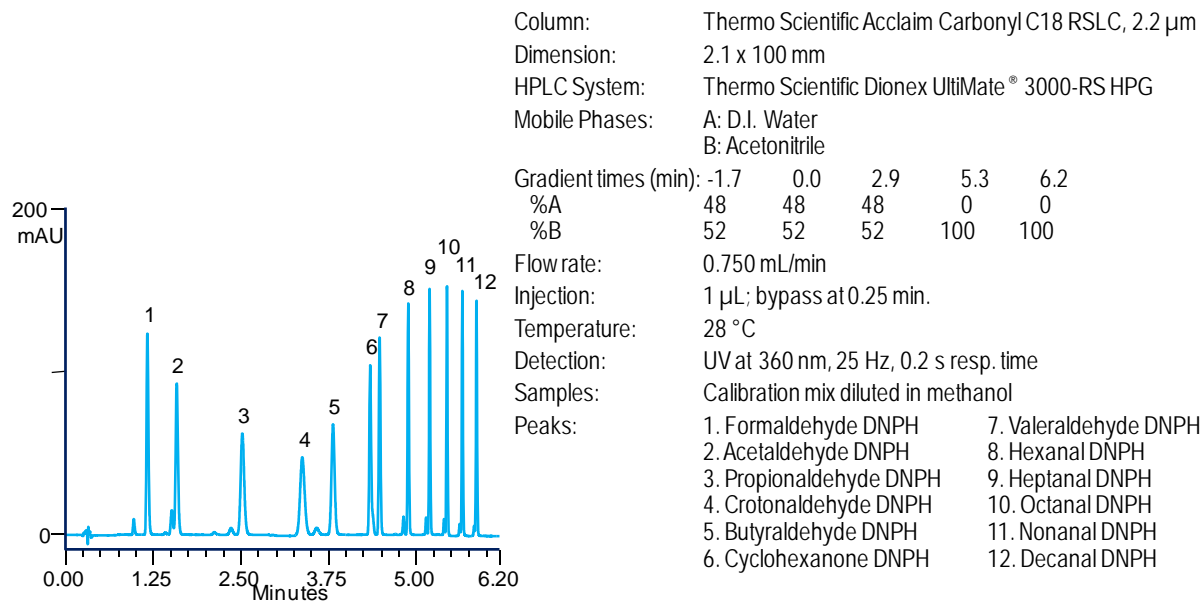


**Figure 3** Rapid CARB-1004 DNPH Aldehyde and Ketone Standards using Acclaim Carbonyl RSLC with Acetonitrile Mobile Phase

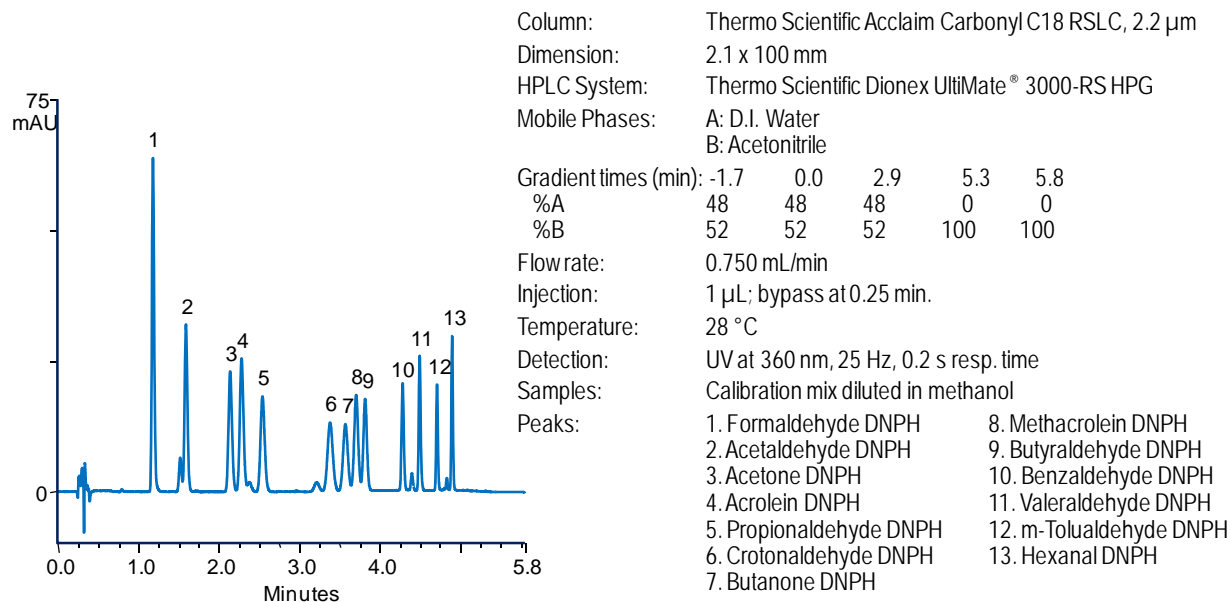


**Figure 4** Rapid EPA 8315 DNPH Aldehyde and Ketone Standards using Acclaim Carbonyl RSLC with Acetonitrile Mobile Phase

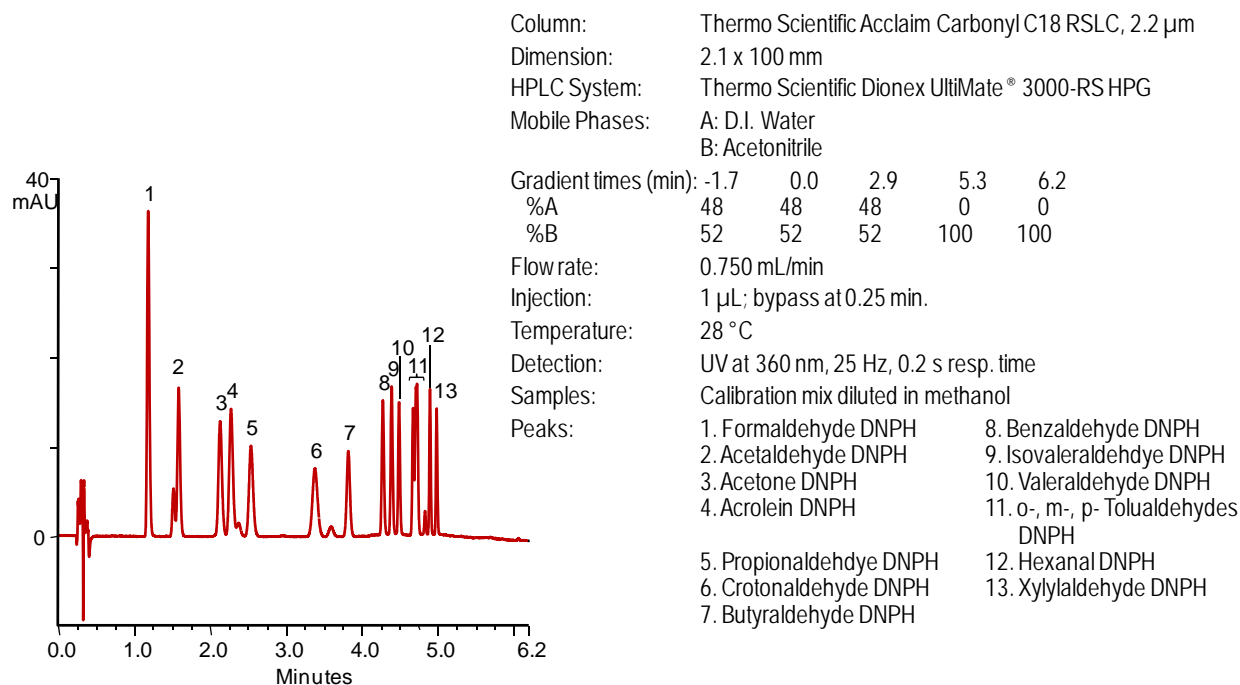


**Figure 5 EPA 1667 DNPH Aldehyde and Ketone Standards using Acclaim Carbonyl RSLC with Acetonitrile Mobile Phase**

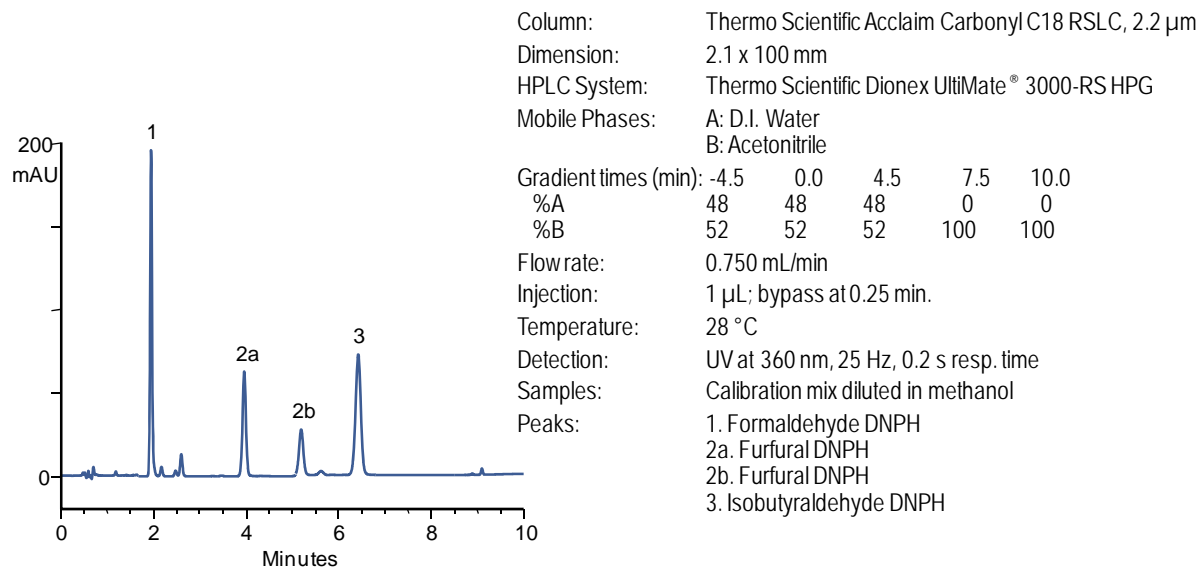


**Figure 6 High-Resolution DNPH Aldehyde and Ketone Standards using Acclaim Carbonyl RSLC with Acetonitrile Mobile Phase**

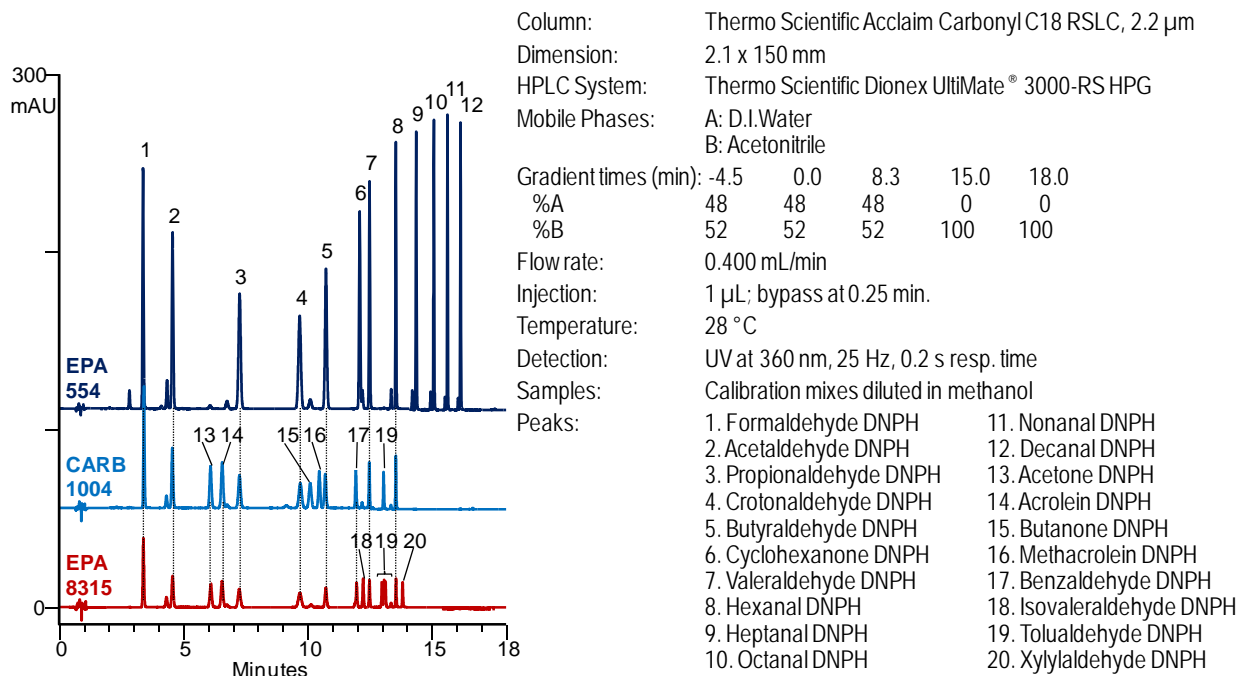

**Figure 7 EPA 554 Carbonyl-DNPH Standards using Acclaim Carbonyl RSLC with Methanol Mobile Phase**

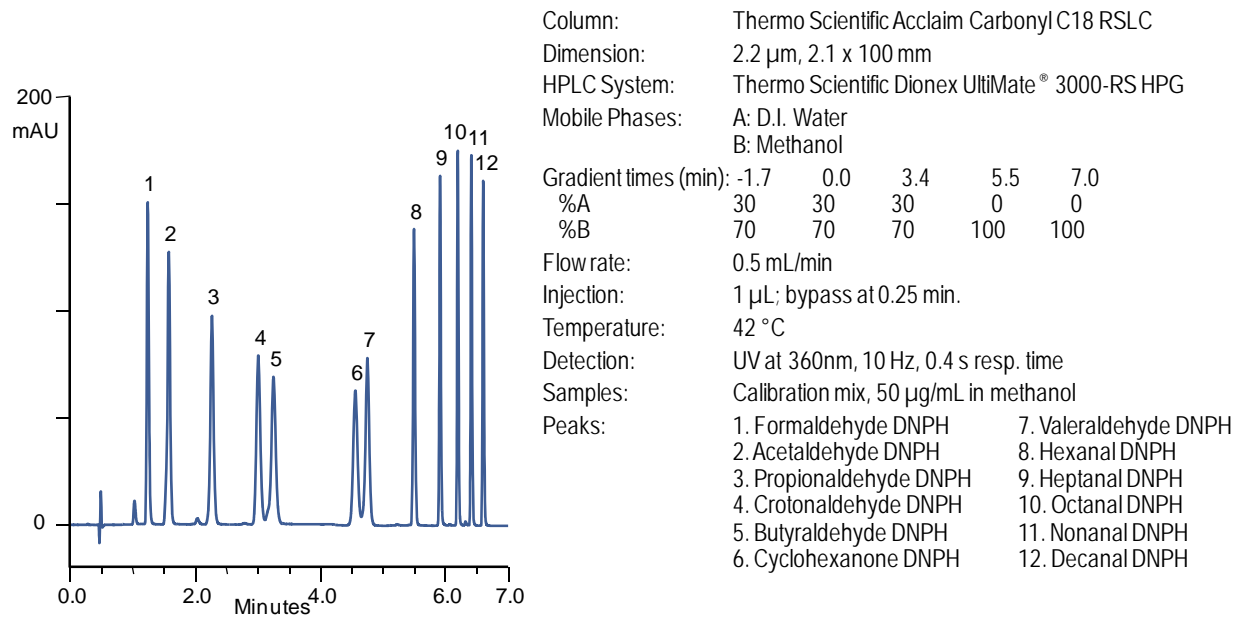

**Figure 8** DNPH Aldehyde and Ketone Standards using Acclaim Carbonyl RSLC with Acetonitrile Mobile Phase

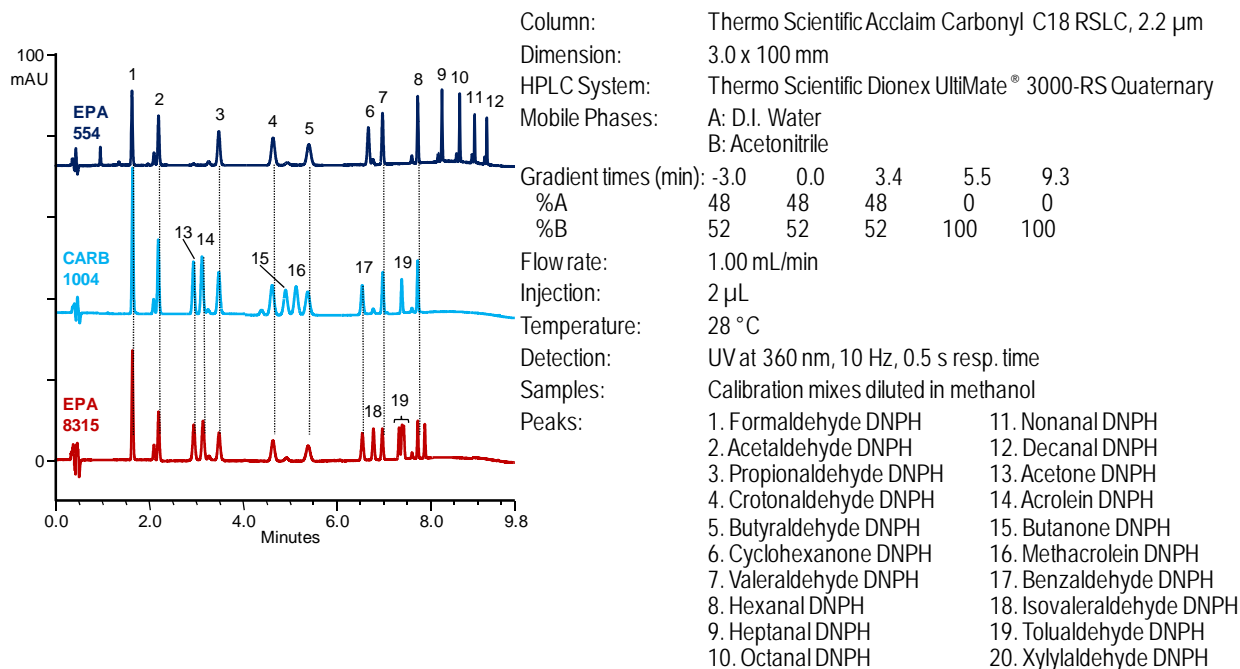



**Figure 9 EPA 554 Carbonyl-DNPH Standards using Acclaim Carbonyl RSLC with Methanol Mobile Phase**


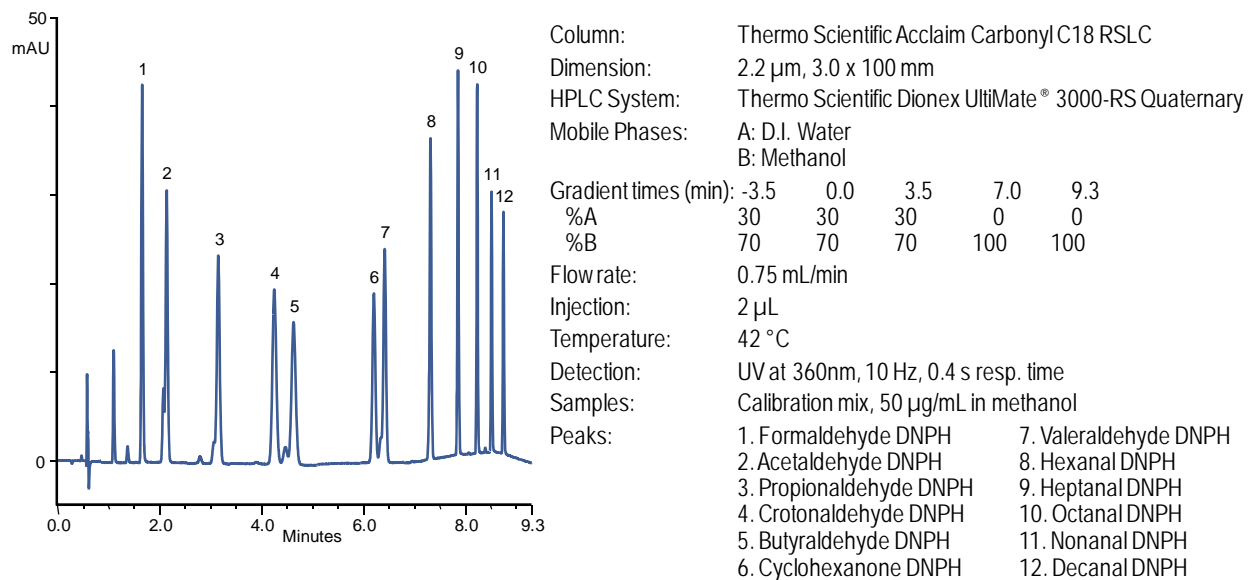


Figure 10 CARB1004 Carbonyl-DNPH Standards using Acclaim Carbonyl 5  $\mu$ m, 4.6 x 250 mm

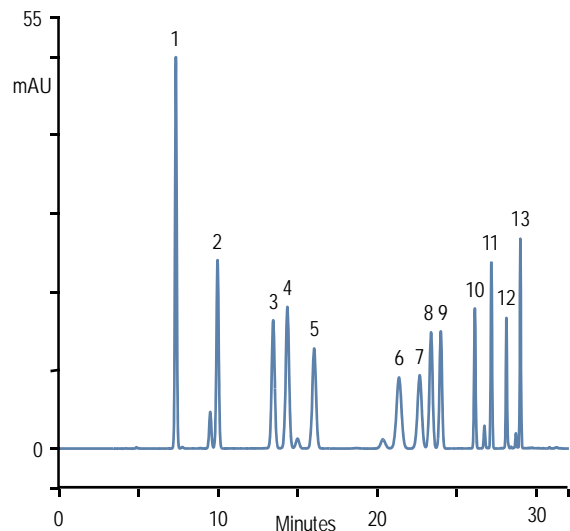

## HPLC Conditions:

Column: Thermo Scientific Acclaim Carbonyl C18 RSLC  
5 $\mu$ m 4.6x250mm

Pump: Thermo Scientific Dionex UltiMate LPG-3400A

Mobile Phases:  
A: Acetonitrile  
B: D.I. Water

Gradient:	-9	0	20	28	32
%A	52	52	52	90	90
%B	48	48	48	10	10

Flowrate: 1.5 mL/min

Temperature: TCC-3100 thermostat at 28 °C

Injection: WPS-3000 SL sampler, 10  $\mu$ L


Detection: VWD-3400, UV 360 nm

Sample: Standards 1.5  $\mu$ g/mL in acetonitrile

## Peaks:

- |                         |                         |
|-------------------------|-------------------------|
| 1. Formaldehyde DNPH    | 8. Methacrolein DNPH    |
| 2. Acetaldehyde DNPH    | 9. Butyraldehyde DNPH   |
| 3. Acetone DNPH         | 10. Benzaldehyde DNPH   |
| 4. Acrolein DNPH        | 11. Valeraldehyde DNPH  |
| 5. Propionaldehyde DNPH | 12. m-Tolualdehyde DNPH |
| 6. Crotonaldehyde DNPH  | 13. Hexaldehyde DNPH    |
| 7. Butanone DNPH        |                         |

Figure 11 CARB1004 Carbonyl-DNPH Standards using Acclaim Carbonyl 3  $\mu$ m, 3 x 150 mm

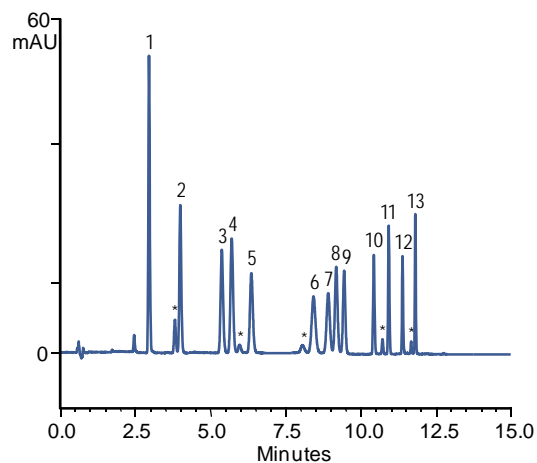

## HPLC Conditions

Column: Thermo Scientific Acclaim Carbonyl C18 RSLC  
3  $\mu$ m, 3 x 150mm

Mobile Phases: A: Acetonitrile  
B: D.I. Water

Pump: Thermo Scientific Dionex UltiMate HPG-3400RS

Gradient:	-3.5	0.0	7.5	11.5	16.0
%A	53	53	53	90	90
%B	47	47	47	10	10

Flow rate: 0.90 mL/min

Pressure: 310 bar max.

Temperature: 28 °C

Injection: 5  $\mu$ L

Detection: UV 360 nm, 10 Hz, 0.5 sec.

Sample: Standards in acetonitrile

## Peaks:

- |                         |                         |
|-------------------------|-------------------------|
| 1. Formaldehyde DNPH    | 8. Methacrolein DNPH    |
| 2. Acetaldehyde DNPH    | 9. Butyraldehyde DNPH   |
| 3. Acetone DNPH         | 10. Benzaldehyde DNPH   |
| 4. Acrolein DNPH        | 11. Valeraldehyde DNPH  |
| 5. Propionaldehyde DNPH | 12. m-Tolualdehyde DNPH |
| 6. Crotonaldehyde DNPH  | 13. Hexaldehyde DNPH    |
| 7. Butanone DNPH        | *. Isomers              |

Figure 12 EPA 8315 Carbonyl-DNPH Standards using Acclaim Carbonyl 3  $\mu$ m, 3 x 150 mm

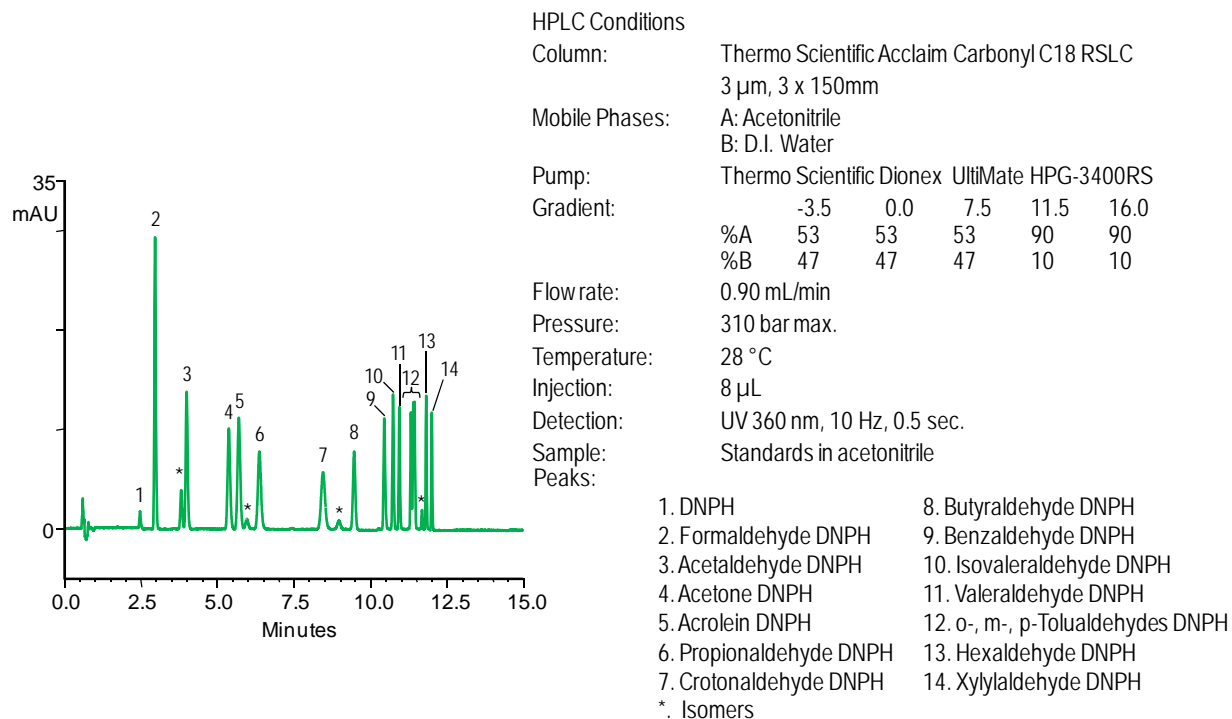


Figure 13 EPA 554 Carbonyl-DNPH Standards using Acclaim Carbonyl 3  $\mu$ m, 3 x 150 mm with Acetonitrile Mobile Phase


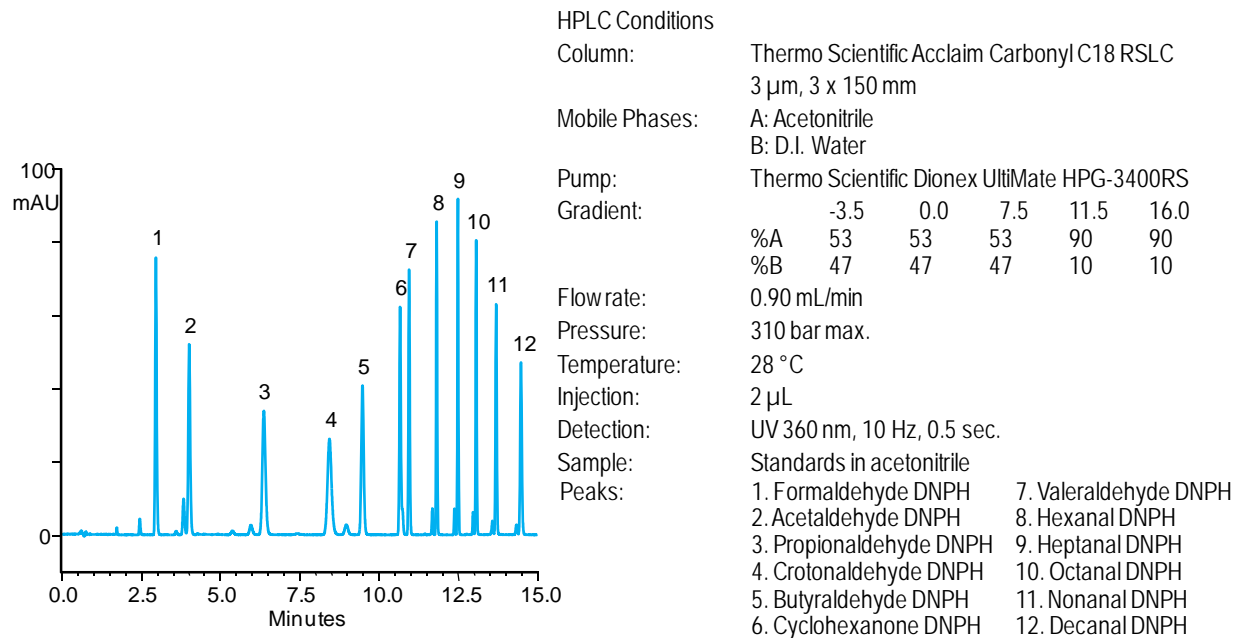


Figure 14 EPA 554 Carbonyl-DNPH Standards using Acclaim Carbonyl 3  $\mu\text{m}$ , 3 x 150 mm with Methanol Mobile Phase

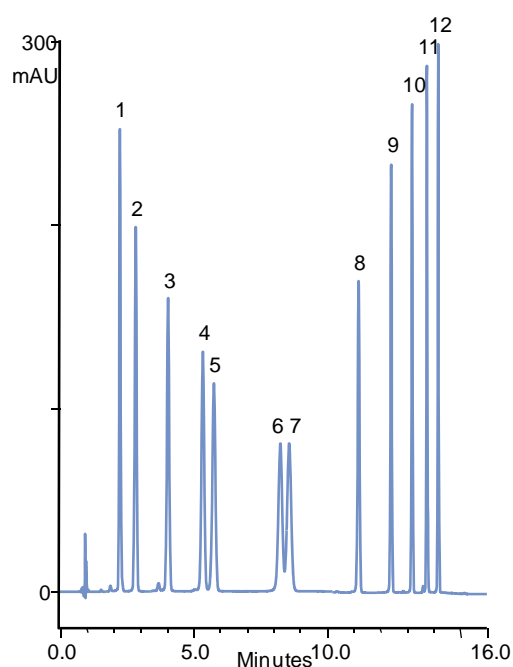

Column: Thermo Scientific Acclaim Carbonyl C18 RSLC  
 Dimension: 3.0  $\mu\text{m}$ , 3.0 x 150 mm  
 HPLC System: UltiMate<sup>®</sup> 3000-RS  
 Mobile Phases: A: D.I Water  
 B: Methanol


Gradient times (min):	-6.0	0.0	7.0	12.5	16.0
%A	28	28	28	0	0
%B	72	72	72	100	100

Flowrate: 0.8 mL/min  
 Injection: 5  $\mu\text{L}$   
 Temperature: 40  $^{\circ}\text{C}$   
 Detection: UV at 230nm, 5 Hz, 0.4 s resp. time  
 Samples: Calibration mix, 50  $\mu\text{g}/\text{mL}$  in ethanol

Peaks:

1. Formaldehyde DNPH	7. Valeraldehyde DNPH
2. Acetaldehyde DNPH	8. Hexanal DNPH
3. Propionaldehyde DNPH	9. Heptanal DNPH
4. Crotonaldehyde DNPH	10. Octanal DNPH
5. Butyraldehyde DNPH	11. Nonanal DNPH
6. Cyclohexanone DNPH	12. Decanal DNPH

Figure 15 EPA 1667 Aldehyde-DNPH Standards using Acclaim Carbonyl 3  $\mu\text{m}$ , 3 x 150 mm

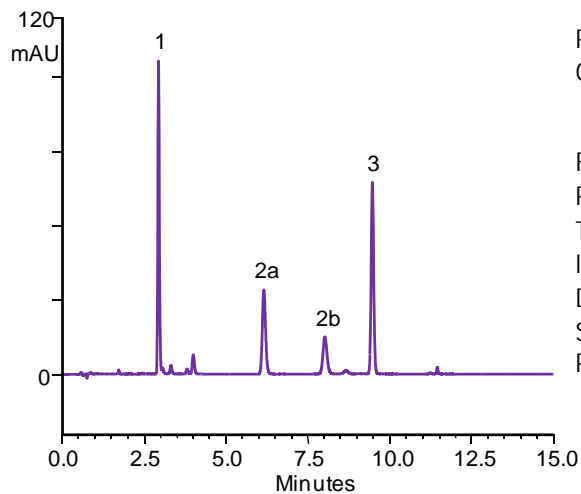

## HPLC Conditions

Column: Thermo Scientific Acclaim Carbonyl C18 RSLC  
3  $\mu\text{m}$ , 3 x 150mm

Mobile Phases: A: Acetonitrile  
B: D.I. Water

Pump: Thermo Scientific Dionex UltiMate HPG-3400RS

Gradient:	-3.5	0.0	7.5	11.5	16.0
%A	53	53	53	90	90
%B	47	47	47	10	10

Flow rate: 0.90 mL/min

Pressure: 310 bar max.

Temperature: 28 °C

Injection: 2  $\mu\text{L}$

Detection: UV 360 nm, 10 Hz, 0.5 sec.

Sample: Standards in acetonitrile

Peaks:  
1. Formaldehyde DNPH  
2a. Furfural DNPH  
2b. Furfural DNPH  
3. Isobutyraldehyde DNPH