



RNA isolation

User manual

NucleoSpin® RNA XS

January 2020 / Rev. 10



www.mn-net.com



RNA isolation

Protocol at a glance (Rev. 10)

XS

NucleoSpin® RNA XS

		I			
1	Supply sample			up to 10 ⁵ cultured cells or 5 mg tissue samples	
2	Lyse and homogenize cells			100 μL RA1 2 μL TCEP	
		V		Mix	
3	Add Carrier RN		5 μL Carrier RNA working solution		
		V	Mix		
4	Filtrate lysate (optional)			11,000 x <i>g</i> , 30 s	
5	Adjust RNA		100 μL 70 % ethanol		
	binding conditions		Mix		
6	Bind RNA		⇒ ►	Load lysate	
				11,000 x <i>g</i> , 30 s	
7	Desalt silica membrane	8	.	100 μL MDB	
				11,000 x <i>g</i> , 30 s	
8	Digest DNA		25 μL DNase reaction mixture		
			RT, 15 min		
9	Wash and dry	8	1st wash	100 μL RA2	
	silica membrane		2 nd wash	400 μL RA3	
			3 rd wash	200 μL RA3	
		1st and 2nd		11,000 x <i>g</i> , 30 s	
		3 rd		11,000 x <i>g</i> , 2 min	
9	Elute highly pure RNA		<u> </u>	10 μL RNase-free H ₂ O	
				11,000 x <i>g</i> , 30 s	



RNA isolation

Table of contents

1	Com	ponents	4
	1.1	Kit contents	4
	1.2	Reagents, consumables, and equipment to be supplied by user	5
	1.3	About this user manual	5
2	Prod	uct description	6
	2.1	The basic principle	6
	2.2	Kit specifications	6
	2.3	Handling, preparation, and storage of starting materials	9
	2.4	Elution procedures	10
	2.5	Stability of isolated RNA	10
3	Stora	age conditions and preparation of working solutions	11
4	Safe	ty instructions	13
5	Proto	ocols	15
	5.1	RNA purification from cultured cells, laser captured cells, or microdissected cryosections	15
	5.2	RNA purification from tissue	18
	5.3	Clean up and concentration of RNA	21
	5.4	rDNase digestion in the eluate	23
6	Appe	endix	25
	6.1	Troubleshooting	25
	6.2	Ordering information	29
	6.3	References	30
	6.4	Product use restriction / warranty	31

1 Components

1.1 Kit contents

	Nucleo Spin® BNA VS			
	NucleoSpin [®] RNA XS			
REF	10 preps 740902.10	50 preps 740902.50	250 preps 740902.250	
Lysis Buffer RA1	6 mL	25 mL	125 mL	
Wash Buffer RA2	2 x 1 mL	15 mL	2 x 15 mL	
Wash Buffer RA3 (Concentrate)*	6 mL	12 mL	50 mL	
Membrane Desalting Buffer MDB	10 mL	10 mL	50 mL	
Reaction Buffer for rDNase	7 mL	7 mL	30 mL	
rDNase, RNase-free (lyphilized)*	1 vial (size A)	1 vial (size C)	2 vials (size D)	
Carrier RNA*	300 μg	300 μg	300 μg	
Reducing Agent TCEP*	14 mg	3 x 14 mg	2 x 107 mg	
RNase-free H ₂ O	13 mL	13 mL	13 mL	
NucleoSpin® Filters (violet rings)	10	50	250	
NucleoSpin [®] RNA XS Columns (light blue rings – plus Collection Tubes)	10	50	250	
Collection Tubes (2 mL)	30	150	750	
Collection Tubes (1.5 mL)	10	50	250	
User manual	1	1	1	

^{*} For preparation of working solutions and storage conditions see section 3.

1.2 Reagents, consumables, and equipment to be supplied by user

Reagents

- 96–100 % ethanol (to prepare Wash Buffer RA3 and for the clean up procedure, section 5.3)
- 70 % ethanol (to adjust RNA binding condition)

Consumables

- 1.5 mL microcentrifuge tubes
- Sterile RNase-free tips

Equipment

- Manual pipettors
- · Centrifuge for microcentrifuge tubes
- Vortex mixer
- Personal protection equipment (lab coat, gloves, goggles)

1.3 About this user manual

It is strongly recommended reading the detailed protocol sections of this user manual if the **NucleoSpin® RNA XS** kit is used for the first time. Experienced users, however, may refer to the Protocol at a glance instead. The Protocol at a glance is designed to be used only as a supplemental tool for quick referencing while performing the purification procedure.

All technical literature is available on the internet at www.mn-net.com.

Please contact Technical Service regarding information about changes of the current user manual compared to previous revisions.

2 Product description

2.1 The basic principle

One of the most important aspects isolating RNA is to prevent degradation of the RNA during the isolation procedure. With the **NucleoSpin® RNA** methods, cells are lysed by incubation in a solution containing large amounts of chaotropic ions. This lysis buffer immediately inactivates RNases – which are present in virtually all biological materials – and creates appropriate binding conditions which favor adsorption of RNA to the silica membrane. Contaminating DNA, which is also bound to the silica membrane, is removed by an rDNase solution which is directly applied onto the silica membrane during the preparation (RNase-free rDNase is supplied with the kit). Simple washing steps with two different buffers remove salts, metabolites and macromolecular cellular components. Pure RNA is finally eluted under low ionic strength conditions with RNase-free H₂O (supplied).

The RNA preparation using **NucleoSpin® RNA** kits can be performed at room temperature. The eluate, however, should be treated with care because RNA is very sensitive to trace contaminations of RNases, often found on general lab ware, fingerprints and dust. To ensure RNA stability, keep RNA frozen at -20 °C for short term or -70 °C for long term storage.

2.2 Kit specifications

- The NucleoSpin® RNA XS kit is recommended for the isolation of RNA from very small samples. Typical sample material comprises small amounts of cells (up to 1 x 10⁵) and tissue (up to 5 mg) such as pellets of cultured cells, laser-captured cells, microdissected cryosections, biopsy samples, fine needle aspirates, and flow cytometer sorted cells (Table 1, page 8).
- The innovative column design with a funnel shaped thrust ring and a small silica membrane area allows elution of RNA in as little as 5–30 μL. Thus, highly concentrated RNA is eluted, ready for common downstream applications (e.g., RT-PCR).
- The RNA yield strongly depends on the sample type, quality, and amount (see Table 2, page 8 for details).
- High quality RNA (RNA Integrity Number (RIN) > 9 according to Agilent 2100
 Bioanalyzer assays) can be obtained from small samples (e.g., 10³ cells, 0.1 mg tissue)
 as well as from larger samples (10⁵ cells, 5 mg tissue). rRNA ratios (28S / 18S) of
 1.8–2.0 can be obtained. Since RNA quality always depends on the sample quality, see
 section 6.3 for further aspects.
- The NucleoSpin® RNA XS kit allows purification of RNA with an A₂₆₀/A₂₈₀ ratio generally exceeding 1.9 (measured in TE buffer pH 7.5). Due to the high RNA purity large amounts of eluates can be used as template in RT-PCR without inhibition (e.g., 8 μL of 10 μL eluates as template in a 20 μL qRT-PCR setup generating stronger signal

compared to reactions with less template in a LightCycler PCR with the Sigma SYBR Green Quantitative RT-PCR Kit).

- The **preparation time** is approximately 45 min for 12 samples.
- As Reducing Agent TCEP (Tris(2-carboxyethyl)phosphine) is supplied in the kit. TCEP is odorless, more stable, more specific for disulfide-bonds, and less toxic than other commonly used reducing agents.
- Carrier RNA (poly(-A) RNA: poly(A) potassium salt, prepared from ADP with polynucleotide phosphorylase) is included for optimal performance with smallest samples.

It is recommended adding Carrier RNA to the sample lysate (20 ng per sample). Such small amounts typically do not interfere with subsequent RT-PCR, even in oligo-dT primed reverse transcriptions. The small amount of Carrier RNA transfered into a reverse transcription reaction is commonly not significantly influencing the outcome of the reaction, due to the large excess of oligo-dT primer. The benefit of adding Carrier RNA to the sample lysate depends on sample type, amount, and kind of downstream RNA analysis. Carrier RNA should be omitted when subsequent to RNA isolation

- a poly-A RNA isolation is performed
- RNA sequencing is performed.
- rDNase is supplied in the kit. DNA contaminations are removed by on column digestion
 with rDNase. For most demanding applications (e.g., expression analysis of plasmid
 transfected cells, plastidial or mitochondrial genes) a subsequent digestion with rDNase
 in the eluate is possible.

Table 1: Kit specifications at a glance			
Parameter	NucleoSpin® RNA XS		
Format	Mini spin column – XS design		
Sample material	Small amounts of tissue < 5 mg tissue, < 100 000 cultured cells		
Fragment size	> 200 nt		
Typical yield	See Table 2 for examples		
A ₂₆₀ /A ₂₈₀	1.9–2.1		
Typical RIN (RNA integrity number)	> 9 (depending on sample quality)		
Elution volume	5–30 μL		
Preparation time	35 min/6 preps		
Binding capacity	110 µg		

Table 2: Overview on average yields of RNA isolation using NucleoSpin® RNA XS				
Sample	Average yield			
10 ⁵ HeLa cells	1000–1500 ng			
10 ⁴ HeLa cells	100–150 ng			
10 ³ HeLa cells	10–15 ng			
10 ² HeLa cells	0.1–1.5 ng			
5 mg mouse kidney	5–8 μg			
1 mg mouse kidney	2 µg			

2.3 Handling, preparation, and storage of starting materials

Work environment

Maintain an RNase-free work environment. Wear gloves at all times during the preparation. Change gloves frequently

Sample storage and RNase inhibition

RNases can rapidly degrade RNA within the samples if samples are not protected from RNase activity after harvest. Use one of the following methods to avoid RNA degradation:

- Use freshly harvested sample for immediate lysis and RNA purification.
- Submerge and store samples in NucleoProtect[®] RNA or similar stabilization solutions.
 Make sure to allow for complete permeation of the sample with the stabilization solution before freezing it. Remove excess stabilization solution from the sample prior to RNA isolation according to the stabilization solution user manual.
- Flash freeze sample in liquid N2 immediately upon harvest and store at -70 °C. Frozen samples are stable up to 6 months. Mortar and pestle can be used to pulverize the sample in a frozen state. Make sure that the sample does not thaw prior to contact with lysis buffer.
- Store samples in Lysis Buffer RA1 after disruption at -70 °C for up to one year, at 4 °C for up to 24 hours or at room temperature for up to several hours. Samples frozen in Lysis Buffer RA1 should be thawed slowly before starting with the isolation of RNA.

Disruption and homogenization of sample material

Cultured cells in suspension:

Collect cells by centrifugation, remove supernatant and immediately add Lysis Buffer RA1 according to step 2 of the standard protocol (see sections 5.1, 6.1).

· Adherent cell cultures (lysis in culture dish):

Completely aspirate cell-culture medium. Immediately add Lysis Buffer RA1 to the cell-culture dish. Avoid incomplete removal of the cell-culture medium in order to allow full lysis activity of the lysis buffer. Continue with lysate filtration (step 3 of the standard protocol).

Adherent cell cultures (lysis after trypsinization):

Aspirate cell-culture medium and wash cells once with PBS. Aspirate PBS. Add $0.1-0.3\,\%$ trypsin in PBS and incubate for an appropriate time to detach the cells from the dish surface. After cell detachment, add medium, transfer cells to an appropriate tube (not supplied), and pellet by centrifugation for 5 min at 300 x g. Remove supernatant and continue with the addition of Lysis Buffer RA1 to the cell pellet.

Animal tissues:

It is essential for efficient RNA preparation that all the RNA contained in the sample is released from the cells by disruption and that the viscosity of the sample is reduced by homogenization.

Thawing of undisrupted animal tissue should only be done in the presence of Buffer RA1 under simultaneous mechanical disruption, for example with a rotor-stator homogenizer or a bead mill. This ensures that the RNA is not degraded by RNases before the preparation has started.

Commonly used techniques for disruption of animal tissues are, for example grinding with **pestle and mortar** or using **a syringe and needle** for multiple passage of the sample through the needle. However, due to the small size of samples to be processed with **NucleoSpin® RNA XS** these disruption methods are often not suitable.

2.4 Elution procedures

A high RNA concentration in the elution fraction is desirable for all typical downstream applications. In particular with regard to limited volumes of reaction mixtures, high RNA concentration can be a crucial criterion. Due to a high default elution volume, standard kits often result in weakly concentrated RNA, if only small samples are processed.

Such RNA often even requires a subsequent concentration to be suitable for the desired application.

In contrast to standard kits, **NucleoSpin® RNA XS** allows an efficient elution in a very small volume resulting in highly concentrated RNA.

Elution volumes in the range of 5-30 µL are recommended, the default volume is 10 µL.

2.5 Stability of isolated RNA

Eluted RNA should immediately be put and always kept on ice during work for optimal stability! Contamination with almost omnipresent RNases (general lab ware, fingerprints, dust) may be a risk for isolated RNA. For short term storage freeze at -20 $^{\circ}$ C, for long term storage freeze at -70 $^{\circ}$ C.

^{*} Add TCEP optional before or after freezing.

3 Storage conditions and preparation of working solutions

Attention: Buffers RA1, RA2, and MDB contain chaotropic salt and detergents. Wear gloves and goggles!

CAUTION: Buffers RA1, RA2, and MDB contain guanidinium thiocyanate which can form highly reactive compounds when combined with bleach (sodium hypochlorite). DO NOT add bleach or acidic solutions directly to the sample-preparation waste.

- Store lyophilized rDNase, Reducing Agent TCEP, and Carrier RNA at 4 °C on arrival (stable up to 1 year).
- All other kit components should be stored at room temperature (18–25 °C) and are stable up to one year. Storage at lower temperatures may cause precipitation of salts.
- Check that 70% ethanol is available as additional solution in the lab to adjust RNA binding conditions in the Buffer RA1 lysate.
- Check that 96–100 % ethanol is available (necessary for clean up protocol only).

Before starting with any NucleoSpin® RNA XS protocol prepare the following:

- rDNase: Add indicated volume (see following table or label on the rDNase vial) of RNase-free H₂O to the rDNase vial and incubate for 1 min at room temperature. Gently swirl the vials to completely dissolve the rDNase. Be careful not to mix rDNase vigorously as rDNase is sensitive to mechanical agitation. Dispense into aliquots and store at -20 °C. The frozen working solution is stable for 6 months. Do not freeze/thaw the aliquots more than three times.
- Reducing Agent TCEP: Add indicated volume of RNase-free H₂O to the TCEP vial and incubate for several minutes at room temperature. Mix the vial to completely dissolve the TCEP. Store dissolved TCEP at -20 °C.
- Carrier RNA: Prepare a stock solution before first time using: Dissolve the Carrier RNA* in 750 μL Buffer RA1 to obtain a 400 ng/μL stock solution. Prepare a working solution before RNA extraction: Dilute 1:100 with Buffer RA1 (e.g., 1 μL Carrier RNA stock solution + 99 μL Buffer RA1) to obtain the working solution of 4 ng/μL. Add 5 μL of this working solution (20 ng) to every lysate (protocol step 3 in section 5). Store stock solution at -20 °C; do not store working solution, prepare it freshly immediately before use.
- Wash Buffer RA3: Add the indicated volume of 96–100 % ethanol to Wash Buffer RA3
 Concentrate. Mark the label of the bottle to indicate that ethanol was added. Store
 Wash Buffer RA3 at room temperature (18–25 °C) for up to one year.

^{*} Due to the production procedure, lyophilized Carrier RNA might hardly be visible in the vial.

	NucleoSpin [®] RNA XS			
REF	10 preps 740902.10	50 preps 740902.50	250 preps 740902.250	
Wash Buffer RA3 (Concentrate)	6 mL Add 24 mL ethanol	12 mL Add 48 mL ethanol	50 mL Add 200 mL ethanol to each bottle	
rDNase, RNase-free (lyophilized)	1 vial (size A) Add 55 μL RNase-free H ₂ O	1 vial (size C) Add 230 μ L RNase-free $\rm H_2O$	2 vials (size D) Add 540 μ L RNase-free H_2 O to each vial	
Carrier RNA	300 μg	300 μg	300 μg	
	Add 750 μL Buffer RA1 to obtain concentrated stock solution. Dilute 1:100 with Buffer RA1 to obtain working solution.			
Reducing Agent TCEP	14 mg Add 100 μL RNase-free H ₂ O	3 x 14 mg Add 100 μL RNase-free H₂O to each vial	$2 \times 107 \text{ mg}$ Add 750 μ L RNase-free H ₂ O to each vial	

4 Safety instructions

The following components of the NucleoSpin® RNA XS kits contain hazardous contents.

Wear gloves and goggles and follow the safety instructions given in this section.

GHS classification

Only harmful features do not need to be labeled with H and P phrases up to 125 mL or 125 g. Mindergefährliche Eigenschaften müssen bis 125 mL oder 125 g nicht mit H- und P-Sätzen gekennzeichnet werden.

Component	Hazard contents	GHS symbol	Hazard phrases	Precaution phrases
Inhalt	Gefahrstoff	GHS-Symbol	H-Sätze	P-Sätze
MDB	Ethanol 5–20 % Ethanol 5–20 %	&	226	210
	CAS 593-84-0, 64-17-5d	WARNING ACHTUNG		
RA1	Guanidinium thiocyanate 45–60 % Guanidinthiocyanat 45–60 % CAS 593-84-0	WARNING ACHTUNG	302, 412	264W, 273, 301+312, 330
RA2	Guanidinium thiocyanate 30–45 % + ethanol 20–35 % Guanidinhydrochlorid 30–45 % + Ethanol 20–35 % CAS 593-84-0, 64-17-5	WARNING ACHTUNG	226, 302, 412	210, 264W, 273, 301+312, 330
rDNase	rDNase 90–100 % rDNase 90–100 % CAS 9003-98-9	DANGER GEFAHR	334	261sh, 342+311
TCEP	Tris-(2-carboxyethyl) phosphine hydrochloride, TCEP(•HCl) 70–100 % Tris-(2-carboxyethyl) phosphinhydrochlorid, TCEP(•HCl) 70–100 % CAS 51805-45-9	WARNING ACHTUNG	315, 319	280sh

Hazard phrases

H226 Flammable liquid and vapour.

Flüssigkeit und Dampf entzündbar.

H302 Harmful if swallowed.

Gesundheitsschädlich bei Verschlucken.

H315 Causes skin irritation.

Verursacht Hautreizungen.

H319 Causes serious eye irritation.

Verursacht schwere Augenreizung.

H334 May cause allergy or asthma symptoms or breathing difficulties if inhaled.

Kann bei Einatmen Allergie, asthmaartige Symptome oder Atembeschwerden verursachen.

H412 Harmful to aquatic life with long lasting effects.

Schädlich für Wasserorganismen, mit langfristiger Wirkung.

Precaution phrases

P 210 Keep away from heat, hot surfaces, sparks, open flames and other ignition

sources. No smoking.

Von Hitze, heissen Oberflächen, Funken, offenen Flammen sowie anderen Zündquellenarten

fernhalten. Nicht rauchen.

P261sh Avoid breathing dust/vapors.

Einatmen von Staub/Dampf vermeiden.

P264W Wash with water thoroughly after handling.

Nach Gebrauch mit Wasser gründlich waschen.

P273 Avoid release to the environment.

Freisetzung in die Umwelt vermeiden.

P280sh Wear protective gloves/eye protection.

Schutzhandschuhe/Augenschutz tragen.

P301+312 IF SWALLOWED: Call a POISON CENTER/doctor if you feel unwell.

BEI VERSCHLUCKEN: Bei Unwohlsein GIFTINFORMATIONSZENTRUM/Arzt anrufen.

P330 Rinse mouth.

Mund ausspülen.

P 342+311 If experiencing respiratory symptoms: Call a POISON CENTER/doctor.

Bei Symptomen der Atemwege: GIFTINFORMATIONSZENTRUM/Arzt anrufen.

For further information please see Material Safety Data Sheets (www.mn-net.com). Weiterführende Informationen finden Sie in den Sicherheitsdatenblättern (www.mn-net.com).

The symbol shown on labels refers to further safety information in this section.

Das auf Etiketten dargestellte Symbol weist auf weitere Sicherheitsinformationen dieses Kapitels hin.

5 Protocols

5.1 RNA purification from cultured cells, laser captured cells, or microdissected cryosections

Before starting the preparation:

 Check if TCEP, Carrier RNA, rDNase, and Wash Buffer RA3 were prepared according to section 3.


1 Supply sample

Provide sample such as a pellet of up to 10⁵ cultured cells, laser captured cells or microdissected cryosections in a microcentrifuge tube (not provided).

For appropriate sample amounts see section 2.2.

2 Lyse and homogenize cells

Add 100 μ L Buffer RA1 and 2 μ L TCEP to the cell sample and vortex vigorously (2 x 5 s).

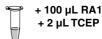

If multiple samples are processed, the preparation of a master-premix is recommended (e.g., 1.1 mL Buffer RA1 and 22 μL TCEP for 10 preparations). Use 102 μL of the premix.

This procedure is usually sufficient to homogenize cultured cells, laser captured cells, or microdissected cryosections. For further comments on homogenization methods see section 2.3.

3 Add Carrier RNA

Add $5 \mu L$ Carrier RNA working solution (20 ng) to the lysate. Mix by vortexing (2 x 5 s). Spin down briefly (approx. 1 s 1000 x g) to clear the lid.



For preparation of Carrier RNA working solution see section 3

4 Filtrate lysate (optional)

Place a **NucleoSpin® Filter** (violet ring) in a Collection Tube (2 mL; supplied), apply the mixture, and centrifuge for **30 s** at **11,000 x** *g*.



11,000 x *g* 30 s

This step may be skipped when working with small amounts of sample, for example less than 10⁵ cells.



Adjust RNA binding condition

Discard the NucleoSpin® Filter (violet ring). Add 100 µL ethanol (70%) to the homogenized lysate and mix by pipetting up and down (5 times).



+ 100 µL 70 % EtOH

Alternatively, add 100 µL ethanol (70%) to the sample in a 1.5 mL microcentrifuge tube (not provided) and mix by vortexing (2 x 5 s). Spin down briefly (approx. 1 s 1000 x a) to clear the lid. Pipette lysate up and down two times before loading the lysate.

Mix

6 **Bind RNA**

For each preparation, take one NucleoSpin® RNA XS Column (light blue ring) placed in a Collection Tube. Load the lysate to the column. Centrifuge for 30 s at 11,000 x q.



Load lysate

Place the column in a new Collection Tube (2 mL).

The maximum loading capacity of NucleoSpin® RNA XS Columns is 600 µL. Repeat the procedure if larger volumes are to be processed.



11,000 x q 30 s

7 Desalt silica membrane

Add 100 uL MDB (Membrane Desalting Buffer) and centrifuge at 11.000 x a for 30 s to dry the membrane. It is not necessary to use a fresh Collection Tube after this centrifugation step.



+ 100 µL MDB

Salt removal will make the following rDNase digest much more effective. If the column outlet has come into contact with the flowthrough for any reason, discard the flowthrough and centrifuge again for 30 s at 11,000 x g.

11,000 x q 30 s

8 **Digest DNA**

reaction mixture in a Prepare **rDNase** microcentrifuge tube (not provided): for each isolation, add 3 μL reconstituted rDNase (also see section 3) to 27 μL Reaction Buffer for rDNase. Mix by flicking the tube.



+ 25 uL rDNase reaction mixture

RT 15 min

Apply 25 uL rDNase reaction mixture directly onto the center of the silica membrane of the column. Close the lid. Incubate at room temperature for 15 min.

It is not necessary to use a new Collection Tube after the incubation step.

9 Wash and dry silica membrane

1st wash

Add 100 μ L Buffer RA2 to the NucleoSpin® RNA XS Column. Incubate for 2 min at RT. Centrifuge for 30 s at 11,000 x g.

Place the column into a new Collection Tube (2 mL).

Buffer RA2 will inactivate the rDNase.

2nd wash

Add **400 µL Buffer RA3** to the NucleoSpin® RNA XS Column. Centrifuge for **30 s** at **11,000 x** *g*. Discard flowthrough and place the column back into the Collection Tube.

Note: Make sure that residual buffer from the previous steps is washed away with Buffer RA3, especially if the lysate has been in contact with the inner rim of the column during loading of the lysate onto the column. For efficient washing of the inner rim flush it with Buffer RA3.

3rd wash

Add **200 µL Buffer RA3** to the NucleoSpin® RNA XS Column. Centrifuge for **2 min** at **11,000 x** *g* to dry the membrane. Place the column into a nuclease-free Collection Tube (1.5 mL; supplied).

If for any reason the liquid level in the Collection Tube has reached the NucleoSpin® RNA XS Column after centrifugation, discard flowthrough and centrifuge again.

+ 100 µL RA2

RT 2 min

11,000 x *g* 30 s



+ 400 µL RA3



11,000 x *g* 30 s

+ 200 µL RA3

11,000 x g 2 min

10 Elute highly pure RNA

Elute the RNA in 10 μ L H_2O (RNase-free; supplied) and centrifuge at 11,000 x g for 30 s.

If higher RNA concentrations or higher elution volumes are desired, elution volume may be varied in the range of 5–30 μ L.

For further details on alternative elution procedures see section 2.4.



_e+ 10 μL RNasefree H₂O

11,000 x *g* 30 s

5.2 RNA purification from tissue

Before starting the preparation:

 Check if TCEP, Carrier RNA, rDNase, and Wash Buffer RA3 were prepared according to section 3.


1 Supply sample

Provide tissue sample such as a biopsy in a microcentrifuge tube (not provided).

For appropriate sample amounts see section 2.2.

2 Lyse and homogenize tissue

Add 200 μ L Buffer RA1 and 4 μ L TCEP to the tissue sample and vortex vigorously (2 x 5 s).

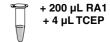

Disruption with a rotor-stator homogenizer or with a shaker and steel balls are recommended methods for the homogenization of tissue samples. For further comments on homogenization methods see section 2.3.



If multiple samples are processed, the preparation of a master-premix is recommended (e.g., 2.2 mL Buffer RA1 and 44 μ L TCEP for 10 preparations). Use 204 μ L of the premix.

3 Add Carrier RNA

Add $5\,\mu$ L Carrier RNA working solution (20 ng) to the lysate. Mix by vortexing (2 x 5 s). Spin down briefly (approx. 1 s 1000 x g) to clear the lid.



+ 5 μL Carrier RNA Mix

For preparation of Carrier RNA working solution see section 3.

4 Filtrate lysate

Reduce viscosity and clear the lysate by filtration through **NucleoSpin® Filter** (violet ring): Place the NucleoSpin® Filter (violet ring) in a Collection Tube (2 mL; provided), apply the mixture, and centrifuge for **30 s** at **11,000 x** *g*.



11,000 x *g* 30 s

In case of visible pellet formation (depending on sample amount and nature), transfer supernatant without any formed pellet to a new 1.5 mL microcentrifuge tube (not included).

5 Adjust RNA binding condition

Discard the NucleoSpin® Filter (violet ring), add 200 µL ethanol (70%) to the homogenized lysate and mix by pipetting up and down (5 times).



+ 200 µL 70 % EtOH

Alternatively, transfer flowthrough into a new 1.5 mL microcentrifuge tube (not provided), add 200 µL ethanol (70%), and mix by vortexing (2 x 5 s). Spin down briefly (approx. 1 s 1000 x g) to clear the lid. Pipette lysate up and down two times before loading the lysate.

After addition of ethanol a stringy precipitate may become visible which will not affect the RNA isolation. Be sure to disaggregate any precipitate by mixing and load all of the precipitate on the column as described in step 6. Do not centrifuge the ethanolic lysate before loading it onto the column in order to avoid pelleting the precipitate.



6 **Bind RNA**

For each preparation, take one NucleoSpin® RNA XS Column (light blue ring) placed in a Collection Tube and load the lysate to the column. Centrifuge for 30 s at 11,000 x q. Place the column in a new Collection Tube (2 mL).



Load lysate

11,000 x q 30 s

The maximum loading capacity of NucleoSpin® RNA XS Columns is 600 μL. Repeat the procedure if larger volumes are to be processed.

7 Desalt silica membrane

Add 100 µL MDB (Membrane Desalting Buffer) and centrifuge at 11,000 x q for 30 s to dry the membrane. It is not necessary to use a fresh Collection Tube after this centrifugation step.



+ 100 µL MDB

11.000 x a 30 s

Salt removal will make the following rDNase digest much more effective. If the column outlet has come into contact with the flowthrough for any reason, discard the flowthrough and centrifuge again for 30 s at 11.000 x g.

8 Digest DNA

Prepare rDNase reaction mixture in a sterile microcentrifuge tube (not provided): for each isolation, add 3 µL reconstituted rDNase (also see section 3) to 27 µL Reaction Buffer for rDNase. Mix by flicking the tube.

+ 25 µL rDNase reaction mixture

> RT 15 min

Apply 25 μ L rDNase reaction mixture directly onto the center of the silica membrane of the column. Close the lid. Incubate at room temperature for 15 min.

It is not necessary to use a new Collection Tube after the incubation step.

9 Wash and dry silica membrane

Add 100 μ L Buffer RA2 to the NucleoSpin® RNA XS Column. Incubate for 2 min at RT. Centrifuge for 30 s at 11,000 x g.

Place the column into a new Collection Tube (2 mL).

Buffer RA2 will inactivate the rDNase.

+ 100 μL RA2 RT 2 min

11,00

11,000 x *g* 30 s

Add $400 \, \mu L$ Buffer RA3 to the NucleoSpin® RNA XS Column. Centrifuge for $30 \, s$ at $11,000 \, x \, g$. Discard flowthrough and place the column back into the Collection Tube.

+ 400 µL RA3

11,000 x *g* 30 s

Add **200 µL Buffer RA3** to the NucleoSpin® RNA XS Column. Centrifuge for **2 min** at **11,000 x** *g* to dry the membrane. Place the column into a nuclease-free Collection Tube (1.5 mL; supplied).

If for any reason the liquid level in the Collection Tube has reached the NucleoSpin® RNA XS Column after centrifugation, discard flowthrough and centrifuge again.



+ 200 µL RA3

2 min

10 Elute highly pure RNA

Elute the RNA in 10 μ L $\rm H_2O$ (RNase-free; supplied) and centrifuge at 11,000 x g for 30 s.

If higher RNA concentrations or higher elution volumes are desired, elution volume may be varied in the range of 5 –30 uL.

For further details on alternative elution procedures see section 2.4



+ 10 μL RNasefree H₂O

> 11,000 x *g* 30 s

5.3 Clean up and concentration of RNA

Before starting the preparation:

Check if Wash Buffer RA3 were prepared according to section 3.

1 Supply sample

Provide **up to 300 µL sample** such as prepurified RNA (e.g., phenol purified) or RNA from reaction mixtures (e.g., labelling reactions) in a microcentrifuge tube (not provided).



Sample

For appropriate sample amounts see section 2.2.

2 Prepare lysis-binding buffer premix

For every 100 μ L of sample combine 25 μ L Buffer RA1 with 75 μ L ethanol (96–100 %) and mix.

If processing multiple samples, the preparation of a master-premix (1 volume Buffer RA1 plus 3 volumes ethanol 96–100%) is recommended.

+ 25 μL RA1 + 75 μL EtOH (96–100%) per 100 μL sample

Mix

3 Add Carrier RNA

Not necessary!

4 Filtrate lysate

Not necessary!

5 Adjust RNA binding condition

Add **one volume of premix** to the sample (e.g., 100 μ L premix to a 100 μ L sample) and mix (2 x 5 s). If necessary, spin down briefly (approx. 1 s 1000 x g) to clear the lid.



Add 1 vol. premix to sample

Mix

6 Bind RNA

For each preparation, take one NucleoSpin® RNA XS Column (light blue ring) placed in a Collection Tube and load the lysate to the column. Centrifuge for 30 s at $11,000 \times g$.



Load lysate

11,000 x *g*

For samples $> 300 \mu L$, load in two steps.

Place the column in a new Collection Tube (2 mL).

For high demanding applications, the recovery rate can be increased as follows: Centrifuge 30 s at 2,000 x g prior to centrifugation for 30 s at 11,000 x g.

7 Desalt silica membrane

Not necessary!

8 Digest DNA

Not necessary!

9 Wash and dry silica membrane

1st wash

Add **400 µL Buffer RA3** to the NucleoSpin® RNA XS Column. Centrifuge for **30 s** at **11,000 x** *g*. Discard flowthrough and place the column back into the Collection Tube.

Ø5

+ 200 µL RA3

+ 400 µL RA3

11,000 x g

30 s

11,000 x *g* 2 min

2nd wash

Add **200 \muL Buffer RA3** to the NucleoSpin® RNA XS Column. Centrifuge for **2 min** at **11,000 x** g to dry the membrane. Place the column into a nuclease-free Collection Tube (1.5 mL; supplied).

If for any reason the liquid level in the Collection Tube has reached the NucleoSpin® RNA XS Column after centrifugation, discard flowthrough and centrifuge again.

10 Elute highly pure RNA

Elute the RNA in 10 μL H_2O (RNase-free; supplied) and centrifuge at 11,000 x g for 30 s.

If higher RNA concentrations or higher elution volumes are desired, elution volume may be varied in the range of 5–30 μ L.



For further details on alternative elution procedures see section 2.4.

5.4 rDNase digestion in the eluate

The on column rDNase digestion in the standard protocol is very efficient and thus results in minimal residual DNA. This DNA will not be detectable in most downstream applications. However, removal of DNA to a completely undetectable level is challenging and the efficiency of an on column DNA digestion is sometimes not sufficient for downstream applications requiring lowest residual content of DNA.

A typical example for such a demanding application is an RT-PCR reaction in which the primer molecules do not differentiate between cDNA (derived from RNA) and contaminating genomic DNA. Especially, if

- high copy number targets are analyzed (e.g., multi gene family, mitochondrial, plastidal or plasmid targets (from transfections))
- the target gene is of a very low expression level
- the amplicon is relatively small (< 200 bp).

DNA digestion in solution can efficiently destroy contaminating DNA. However, stringent RNase control and subsequent repurification of the RNA (in order to remove buffer, salts, DNase and digested DNA) are usually required.

The high quality, recombinant, RNase-free DNase (rDNase) in the **NucleoSpin® RNA XS** kits facilitates such a digestion in solution in order to remove even traces of contaminating DNA.

A Digest DNA (Reaction setup)

Prepare enzyme-buffer premix: Add 1 μ L rDNase to 10 μ L Reaction Buffer for rDNase.

Add 1/10 volume of enzyme-buffer premix to the eluted RNA (e.g., to 10 μ L RNA add 1 μ L of the premix comprising buffer and enzyme).

B Incubate sample

Incubate for 10 min at 37 °C.

C Repurify RNA

Repurify RNA with a suitable RNA cleanup procedure, for example following section 2.4, by ethanol precipitation or with the NucleoSpin® RNA Clean-up XS kit (see ordering information).

Ethanol precipitation, exemplary

Add 0.1 volume of 3 M sodium acetate, pH 5.2 and 2.5 volumes of $96-100\,\%$ ethanol to one volume of sample. Mix thoroughly.

Incubate several minutes to several hours at -20 °C or 4 °C, respectively.

<u>Note:</u> Choose long incubation times if the sample contains low RNA concentration. Short incubation times are sufficient if the sample contains high RNA concentration.

Centrifuge for 10 min at max. speed.

Wash RNA pellet with 70 % ethanol.

Dry RNA pellet and resuspend RNA in RNase-free H₂O.

6 Appendix

6.1 Troubleshooting

Problem Possible cause and suggestions

RNase contamination

RNA is degraded/ no RNA obtained

 Create an RNase free working environment. Wear gloves during all steps of the procedure. Change gloves frequently. Use of sterile, disposable polypropylene tubes is recommended. Keep tubes closed whenever possible during the preparation. Glassware should be oven-baked for at least 2 hours at 250 °C before use.

Reagents not applied or restored properly

- Reagents not properly restored. Add the indicated volume of RNase-free H₂O to rDNase vial and 96 % ethanol to Buffer RA3 Concentrate and mix. Reconstitute and store lyophilized rDNase according to instructions given in section 3.
- Sample and reagents have not been mixed completely. Always vortex vigorously after each reagent has been added.
- No ethanol has been added after lysis. Binding of RNA to the silica membrane is only effective in the presence of ethanol.

Kit storage

Poor RNA quality or yield

- Reconstitute and store lyophilized rDNase according to instructions given in section 3.
- Store other kit components at room temperature. Storage at low temperatures may cause salt precipitation.
- Keep bottles tightly closed in order to prevent evaporation or contamination.

Ionic strength and pH influence A_{260} absorption as well as ratio A_{260}/A_{280}

- For adsorption measurement, use 5 mM Tris pH 8.5 as diluent.
 Please see also:
- Manchester, K L. 1995. Value of A_{260}/A_{280} ratios for measurement of purity of nucleic acids. Biotechniques 19, 208–209.
- Wilfinger, W W, Mackey, K and Chomczyski, P. 1997. Effect of pH and ionic strength on the spectrophotometric assessment of nucleic acid purity. Biotechniques 22, 474–481.

Problem Possible cause and suggestions Sample material Sample material not stored properly. Whenever possible, use fresh material. If this is not possible, store samples in NucleoProtect® RNA stabilization solution or flash freeze the samples in liquid N₂. Store unprotected samples at -70 °C. Never allow unprotected samples to thaw before addition of Buffer RA1. Perform disruption of samples immediately after addition of Lysis Buffer RA1. Poor RNA quality or yield Sample material not stored properly. Whenever possible, use (continued) fresh material. If this is not possible, flash freeze the samples in liquid N2. Samples should always be kept at -70 °C. Never allow tissues to thaw before addition of Buffer RA1. Perform disruption of samples immediately after addition of Lysis Buffer RA1. Insufficient disruption and/or homogenization of starting material. Ensure thorough sample disruption and use NucleoSpin® Filters for easy homogenization of disrupted starting material. Carry-over of guanidinium thiocyanate Carefully load the lysate to the NucleoSpin® RNA XS Column and try to avoid a contamination of the upper part of the column and Low A₂₆₀/A₂₃₀ the column lid. ratio Make sure that residual Wash Buffer RA2 is washed away with Wash Buffer RA3. This may be done by applying Buffer RA3 to the inner rim of the column. Sample material Clogged Too much starting material used. Overloading may lead to NucleoSpin® decreased overall yield. Reduce amount of sample material or use Column/ larger volume of Buffer RA1. Poor RNA quality Insufficient disruption and/or homogenization of starting material. or yield Ensure thorough sample disruption and use NucleoSpin® Filters for easy homogenization of disrupted starting material. rDNase not active Reconstitute and store lyophilized rDNase according to instructions given in section 3. Contamination DNase solution not properly applied of RNA with Pipette rDNase solution directly onto the center of the silica genomic DNA membrane and close the lid. Too much cell material used Reduce quantity of cells or tissue used.

Problem

Possible cause and suggestions

DNA detection system too sensitive

 The amount of DNA contamination is effectively reduced during the on column digestion with rDNase. Anyhow, it can not be guaranteed that the purified RNA is 100 % free of DNA, therefore in very sensitive applications it might still be possible to detect DNA.

Contamination of RNA with genomic DNA (continued)

The probability of DNA detection with PCR increases with:

- the number of DNA copies per preparation: single copy target < plastidial/mitochondrial target < plasmid transfected into cells
- decreasing of PCR amplicon size.
- Use larger PCR targets (e.g., > 500 bp) or intron spanning primers if possible.
- Use support protocol 5.4 for subsequent rDNase digestion in solution.

Carry-over of ethanol or salt

 Do not let the flowthrough touch the column outlet after the second Buffer RA3 wash. Be sure to centrifuge at the corresponding speed for the respective time in order to remove ethanolic Buffer RA3 completely.

Suboptimal performance of RNA in downstream experiments

- Check if Buffer RA3 has been equilibrated to room temperature before use. Washing at lower temperatures lowers efficiency of salt removal by Buffer RA3.
- Depending on the robustness of the used RT-PCR system, RT-PCR might be inhibited if complete eluates are used as template for RT-PCR. Use less eluate as template.

Store isolated RNA properly

 Eluted RNA should always be kept on ice for optimal stability since trace contaminations of omnipresent RNases (general lab ware, fingerprints, dust) will degrade the isolated RNA. For short term storage freeze at -20 °C, for long term storage freeze at -70 °C.

Problem

Possible cause and suggestions

Silica abrasion from the membrane

Discrepancy between A₂₆₀ quantification values and PCR quantification values • Due to the typically low RNA content in very small samples and the resulting low total amount of isolated RNA, an RNA quantification via A₂₆₀ absorption measurement is often hampered due to the low sensitivity of the absorption measurement. When performing absorption measurements close to the detection limit of the photometer, the measurement may be influenced by minor amounts of silica abrasion. In order to prevent incorrect A₂₆₀-quantification of small RNA amounts centrifuge the eluate for 30 s at > 11.000 x g and take an aliquot for measurement without disturbing any sediment. Alternatively, use a silica abrasion insensitive RNA quantification method (e.g., RiboGreen fluorescent dye).

Measurement not in the range of photometer detection limit

Unexpected A₂₆₀/A₂₈₀ ratio

• In order to obtain a significant A_{260}/A_{280} ratio it is necessary that the initially measured A_{260} and A_{280} values are significantly above the detection limit of the photometer used. An A_{280} value close to the background noise of the photometer will cause unexpected A_{260}/A_{280} ratios.

6.2 Ordering information

Product	REF	Pack of
NucleoSpin® RNA XS	740902.10/.50/.250	10/50/250
NucleoSpin® RNA Clean-up XS	740903.10/.50/.250	10/50/250
NucleoSpin® totalRNA FFPE XS	740969.10/.50/.250	10/50/250
NucleoSpin® RNA	740955.10/.50/.250	10/50/250
NucleoZOL	740404.200	200 mL
NucleoSpin® RNA Blood	740200.10/.50	10/50
NucleoSpin® totalRNA FFPE	740982.10/.50/.250	10/50/250
NucleoSpin® RNA Midi	740962.20	20
NucleoSpin® RNA/Protein	740933.10/.50/.250	10/50/250
NucleoSpin® TriPrep	740966.10/.50/.250	10/50/250
NucleoSpin® RNA Clean-up	740948.10/.50/.250	10/50/250
NucleoSpin® RNA/DNA Buffer Set	740944	Suitable for 100 preps
Buffer RA1	740961 740961.500	60 mL 500 mL
rDNase Set	740963	1 set
Reducing Agent TCEP	740395.107	107 mg
NucleoSpin® Filters	740606	50
Collection Tubes (2 mL)	740600	1000
NucleoProtect® RNA	740400.50/.250/.500	50/250/500 mL

Visit www.mn-net.com for more detailed product information.

6.3 References

Fleige S, Pfaffl MW.: RNA integrity and the effect on the real-time qRT-PCR performance. Mol Aspects Med. 2006 Apr-Jun; 27(2–3):126–39. Epub 2006 Feb 15. Review.

Imbeaud S, Graudens E, Boulanger V, Barlet X, Zaborski P, Eveno E, Mueller O, Schroeder A, Auffray C.: Towards standardization of RNA quality assessment using user-independent classifiers of microcapillary electrophoresis traces. Nucleic Acids Res. 2005 Mar 30;33(6):e56.

Miller CL, Diglisic S, Leister F, Webster M, Yolken RH.: Evaluating RNA status for RT-PCR in extracts of postmortem human brain tissue. Biotechniques. 2004 Apr; 36(4):628–33.

Schoor O, Weinschenk T, Hennenlotter J, Corvin S, Stenzl A, Rammensee HG, Stevanovic S.: Moderate degradation does not preclude microarray analysis of small amounts of RNA. Biotechniques. 2003 Dec; 35(6):1192–6, 1198–201.

6.4 Product use restriction / warranty

NucleoSpin® RNA XS kit components are intended, developed, designed, and sold FOR RESEARCH PURPOSES ONLY, except, however, any other function of the product being expressly described in original MACHEREY-NAGEL product leaflets.

MACHEREY-NAGEL products are intended for GENERAL LABORATORY USE ONLY! MACHEREY-NAGEL products are suited for QUALIFIED PERSONNEL ONLY! MACHEREY-NAGEL products shall in any event only be used wearing adequate PROTECTIVE CLOTHING. For detailed information please refer to the respective Material Safety Data Sheet of the product! MACHEREY-NAGEL products shall exclusively be used in an ADEQUATE TEST ENVIRONMENT. MACHEREY-NAGEL does not assume any responsibility for damages due to improper application of our products in other fields of application. Application on the human body is STRICTLY FORBIDDEN. The respective user is liable for any and all damages resulting from such application.

DNA/RNA/PROTEIN purification products of MACHEREY-NAGEL are suitable for IN-VITRO-USES ONLY!

ONLY MACHEREY-NAGEL products specially labeled as IVD are also suitable for IN-VITRO-diagnostic use. Please pay attention to the package of the product. IN-VITRO-diagnostic products are expressly marked as IVD on the packaging.

IF THERE IS NO IVD SIGN, THE PRODUCT SHALL NOT BE SUITABLE FOR IN-VITRO-DIAGNOSTIC USE!

ALL OTHER PRODUCTS NOT LABELED AS IVD ARE NOT SUITED FOR ANY CLINICAL USE (INCLUDING, BUT NOT LIMITED TO DIAGNOSTIC, THERAPEUTIC AND/OR PROGNOSTIC USE).

No claim or representations is intended for its use to identify any specific organism or for clinical use (included, but not limited to diagnostic, prognostic, therapeutic, or blood banking). It is rather in the responsibility of the user or - in any case of resale of the products - in the responsibility of the reseller to inspect and assure the use of the DNA/RNA/protein purification products of MACHEREY-NAGEL for a well-defined and specific application.

MACHEREY-NAGEL shall only be responsible for the product specifications and the performance range of MN products according to the specifications of in-house quality control, product documentation and marketing material.

This MACHEREY-NAGEL product is shipped with documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. MACHEREY-NAGEL's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Supplementary reference is made to the general business terms and conditions of MACHEREY-NAGEL, which are printed on the price list. Please contact us if you wish to get an extra copy.

There is no warranty for and MACHEREY-NAGEL is not liable for damages or defects arising in shipping and handling (transport insurance for customers excluded), or out of accident or improper or abnormal use of this product; defects in products or components not manufactured by MACHEREY-NAGEL, or damages resulting from such non-MACHEREY-NAGEL components or products.

MACHEREY-NAGEL makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING,

WITHOUT LIMITATION, AS TO THE SUITABILITY, REPRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO MACHEREY-NAGEL PRODUCTS.

In no event shall MACHEREY-NAGEL be liable for claims for any other damages, whether direct, indirect, incidental, compensatory, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure of MACHEREY-NAGEL products to perform in accordance with the stated specifications. This warranty is exclusive and MACHEREY-NAGEL makes no other warranty expressed or implied.

The warranty provided herein and the data, specifications and descriptions of this MACHEREY-NAGEL product appearing in MACHEREY-NAGEL published catalogues and product literature are MACHEREY-NAGEL's sole representations concerning the product and warranty. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agent or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.


Product claims are subject to change. Therefore please contact our Technical Service Team for the most up-to-date information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Applications mentioned in MACHEREY-NAGEL literature are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of any of those applications.

Last updated: 07 / 2010, Rev. 03

Please contact:

MACHEREY-NAGEL GmbH & Co. KG

Tel.: +49 24 21 969-270 tech-bio@mn-net.com

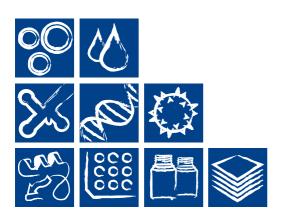

www.mn-net.com



MACHEREY-NAGEL GmbH & Co. KG · Neumann-Neander-Str. 6-8 · 52355 Düren · Germany FR:

DE/International:

CH: Tel.: +49 24 21 969-0 Tel.: +41 62 388 55 00 Fax: +49 24 21 969-199 Fax: +41 62 388 55 05 E-mail: info@mn-net.com

Tel.: +33 388 68 22 68 Fax: +33 388 51 76 88 E-mail: sales-ch@mn-net.com E-mail: sales-fr@mn-net.com

US:

Tel.: +1 484 821 0984 Fax: +1 484 821 1272 E-mail: sales-us@mn-net.com

A034939/0101.5