


MACHEREY-NAGEL

# User manual





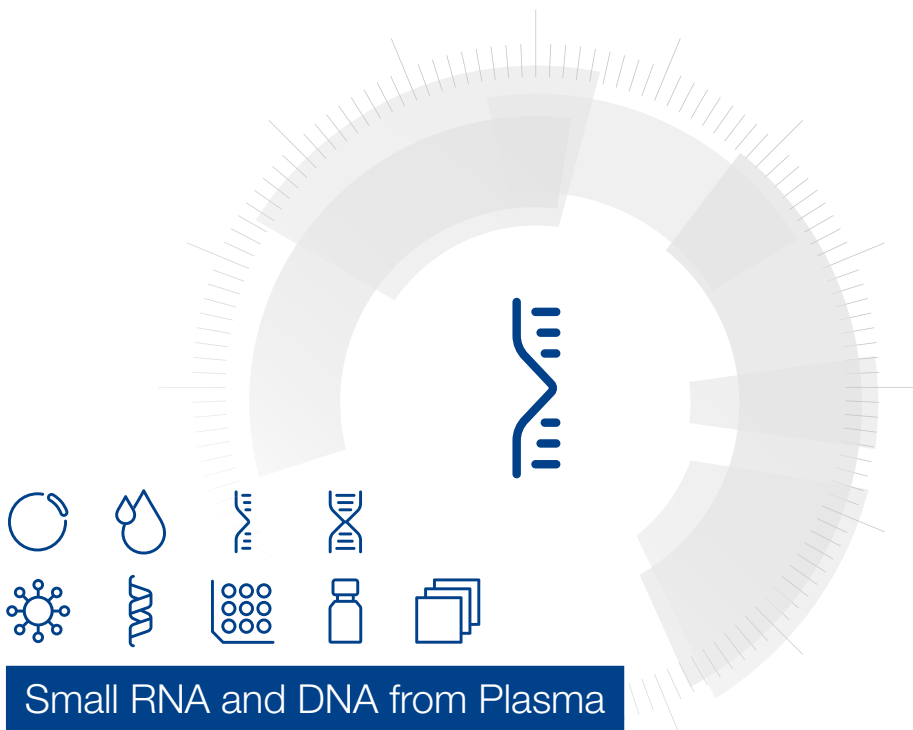











## Small RNA and DNA from Plasma

■ NucleoSpin® miRNA Plasma

April 2022 / Rev. 07
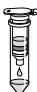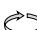
# Isolation of miRNA from plasma

## Protocol at a glance (Rev.07)

<b>NucleoSpin® miRNA Plasma</b>		
<b>1 Prepare sample</b>		300 µL plasma or serum* 90 µL MLP Vortex 5 s RT, 3 min
<b>2 Precipitate protein</b>	 	30 µL MPP Vortex 5 s RT, 1 min 11,000 x g, 3 min
<b>3 Transfer supernatant</b>		Transfer clear supernatant to Collection Tube (2 mL, lid)
<b>4 Adjust binding conditions</b>	 	400 µL isopropanol Vortex 5 s
<b>5 Bind RNA and DNA</b>	 	Load sample on NucleoSpin® miRNA Column RT, 2 min 11,000 x g, 30 s
<b>6 Optional: Digest DNA</b>	 	<b>Optional:</b> <div style="display: flex; flex-direction: column; gap: 10px;"> <div> <span style="background-color: black; color: white; padding: 2px 5px;"><b>1<sup>st</sup></b></span>                     700 µL MW2                      11,000 x g, 30 s                 </div> <div> <span style="background-color: black; color: white; padding: 2px 5px;"><b>2<sup>nd</sup></b></span>                     250 µL MW2                      11,000 x g, 2 min                 </div> </div> 50 µL rDNase in Reaction Buffer for rDNase RT, 15 min
<b>7 Wash and dry</b>	 	<div style="display: flex; flex-direction: column; gap: 10px;"> <div> <span style="background-color: black; color: white; padding: 2px 5px;"><b>1<sup>st</sup></b></span>                     100 µL MW1                      11,000 x g, 30 s                 </div> <div> <span style="background-color: black; color: white; padding: 2px 5px;"><b>2<sup>nd</sup></b></span>                     700 µL MW2                      11,000 x g, 30 s                 </div> <div> <span style="background-color: black; color: white; padding: 2px 5px;"><b>3<sup>rd</sup></b></span>                     250 µL MW2                      11,000 x g, 2 min                 </div> </div>
<b>8 Elute RNA</b>	 	30 µL RNase-free H <sub>2</sub> O RT, 1 min 11,000 x g, 1 min

\* Larger sample volumes can be processed when buffer volumes of MLP, MPP, and isopropanol are increased proportionally. Multiple loading is necessary in step 5 (see section 2.3 for more information).

## Table of contents

1	Components	4
1.1	Kit contents	4
1.2	Reagents, consumables, and equipment to be supplied by user	5
1.3	About this user manual	5
2	Product description	6
2.1	The basic principle	6
2.2	Kit specifications	6
2.3	Amount of starting material	6
2.4	Proteinase K digest	8
2.5	Addition of carrier	8
2.6	Elution procedures	8
3	Storage conditions and preparation of working solutions	9
4	Safety instructions	10
4.1	Disposal	10
5	Protocol - small RNA and DNA purification from plasma or serum	11
6	Appendix	14
6.3	Product use restriction / warranty	15
6.1	Troubleshooting	14
6.2	Ordering information	15
6.3	Product use restriction / warranty	15

# 1 Components

## 1.1 Kit contents

<b>NucleoSpin® miRNA Plasma</b>			
<b>REF</b>	<b>10 preps 740981.10</b>	<b>50 preps 740981.50</b>	<b>250 preps 740981.250</b>
Lysis Buffer MLP	3 mL	13 mL	75 mL
Protein Precipitation Buffer MPP	5 mL	5 mL	25 mL
Reaction Buffer for rDNase	7 mL	7 mL	30 mL
rDNase, RNase-free (lyophilized)	1 vial (size C)	1 vial (size C)	5 vials (size C)
Wash Buffer MW1	10 mL	10 mL	35 mL
Wash Buffer MW2 (Concentrate)*	6 mL	25 mL	100 mL
RNase-free H <sub>2</sub> O	13 mL	13 mL	13 mL
NucleoSpin® miRNA Columns (green rings)	10	50	250
Collection Tubes (1.5 mL)	10	50	250
Collection Tubes (2 mL)	10	50	250
Collection Tubes (2 mL, lid)	10	50	250
User manual	1	1	1

\* For preparation of working solutions and storage conditions see section 3.

## 1.2 Reagents, consumables, and equipment to be supplied by user

### Reagents

- 96 – 100 % ethanol
- Isopropanol

### Consumables

- RNase-free disposable pipette tips

### Equipment

- Manual pipettors
- Vortexer
- Centrifuge for microcentrifuge tubes
- Personal protection equipment (lab coat, gloves, goggles)

## 1.3 About this user manual

It is strongly recommended that first-time users of the **NucleoSpin® miRNA Plasma** kit read the detailed protocol section of this user manual. Experienced users, however, may refer to the Protocol at a glance instead. The Protocol at a glance is designed to be used only as a supplemental tool for quick referencing while performing the purification procedure.

All technical literature is available on the internet at [www.mn-net.com](http://www.mn-net.com).

## 2 Product description

### 2.1 The basic principle

The **NucleoSpin® miRNA Plasma** kit offers the unique feature to isolate total RNA including small RNA and DNA from serum and plasma without the need to resort to the cumbersome phenol/chloroform extraction or a time consuming proteinase digest.

The sample material is denatured in Lysis Buffer MLP. The protein is then precipitated by Protein Precipitation Buffer MPP and pelleted by centrifugation.

After the removal of protein the binding conditions for nucleic acids are adjusted by adding isopropanol.

Total nucleic acids are bound to the NucleoSpin® miRNA Column. Optionally, DNA can be removed by an on-column rDNase digest. The remaining nucleic acids are washed and eluted with minimal amounts of RNase-free water.

### 2.2 Kit specifications

- The **NucleoSpin® miRNA Plasma** kit is designed for the isolation of RNA and DNA from cell free blood plasma or serum.
- **rDNase** is provided for an optional on-column digest to remove traces of genomic DNA.
- The eluted RNA and DNA are ready to use for all standard downstream applications, for example, qPCR, qRT-PCR, Northern Blot, chip hybridization.

**Table 1: Kit specifications at a glance**

Parameter	NucleoSpin® miRNA Plasma
Sample size	300 µL blood plasma or serum
Binding capacity	200 µg
Elution volume	20 – 50 µL
Preparation time	40 min/10 preps (without rDNase digest) 70 min/10 preps (with rDNase digest)

### 2.3 Amount of starting material


The standard procedure allows to process 300 µL of sample material with only one loading step onto the NucleoSpin® miRNA Column. This is usually enough to detect also low abundance miRNA in plasma or serum.

If larger sample volumes are to be used to increase the sensitivity even further, the volumes of Buffer MLP and Buffer MPP as well as the isopropanol have to be increased proportionally. Multiple loading steps per sample are necessary.

But, consider that doubling or tripling the starting volume will result in an only 1 – 1.5 cycles earlier signal in qRT-PCR, which is rather insignificant for detection sensitivity compared to the much larger differences that occur from sample to sample or between different miRNAs.

Furthermore, if plasma quality is low, co-purified RT-PCR inhibitors might require diluting the eluate and thereby counteract the increased yields.

Figure 1 shows qRT-PCR results from eight different plasma samples analyzed for miR-16 (very high expression) and miR-1 (very low expression). They differ in average by about 14 cycles, which is a difference in expression by 3–4 orders of magnitude (1.000 – 10.000 fold)! Doubling the amount of sample material would just shift the mean values from 32.6 to 31.6 (miR-1) and from 18.3 to 17.3 (miR-16). This is much less than the variation from sample to sample and can thus be neglected.

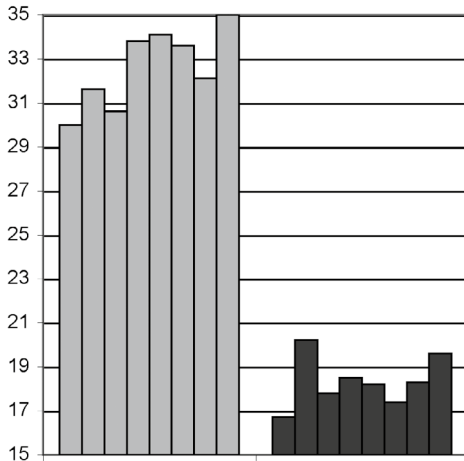

**Figure 1 qRT-PCR quantification of miR-1 and miR-16**

MiRNA from only 300  $\mu\text{L}$  of eight different blood plasma samples were purified and eluted in 30  $\mu\text{L}$  RNase-free water. 2  $\mu\text{L}$  of each eluate were used for a 10  $\mu\text{L}$  RT reaction (Applied Biosystems, TaqMan<sup>®</sup> MicroRNA RT Kit) with miRNA specific primers (Applied Biosystems, hsa-miR-1 MicroRNA Assay and hsa-miR-16 MicroRNA Assay). The RT reaction mix was diluted 1:10. Only 4  $\mu\text{L}$  of the resulting 100  $\mu\text{L}$  were used for the PCR reaction (Applied Biosystems TaqMan<sup>®</sup> Universal PCR Master Mix) in combination with the MicroRNA Assays for specific priming mentioned before.

### 2.3.1 Preparation of plasma from human EDTA blood

- 1 Centrifuge fresh blood sample for 10 min at 2,000 x g.
- 2 Remove the plasma without disturbing sedimented cells.
- 3 Freeze plasma at -20 °C for storage upon RNA isolation.
- 4 Thaw frozen plasma samples prior to RNA isolation and centrifuge for 3 min at  $\geq 11,000$  x g in order to remove residual cells, cell debris, and particulate matter. Use the supernatant for RNA isolation.

## 2.4 Proteinase K digest

A short protein digestion step might increase miRNA yield, especially for low quality, hemolyzed plasma. Add 10  $\mu\text{L}$  of Proteinase K (~ 30  $\mu\text{g}/\mu\text{L}$  in Buffer PB, see ordering information) to 300  $\mu\text{L}$  plasma, incubate for 10 min at 37 °C, and then proceed with addition of Buffer MLP according to the protocol.

## 2.5 Addition of carrier

To improve RNA/DNA binding to the NucleoSpin® miRNA Column, carrier can be added to the sample after the removal of precipitated protein. Slightly higher yields could be found with 2  $\mu\text{g}$  of glycogen or 5  $\mu\text{g}$  of LPA (linear polyacrylamide). However, negligible effects were observed for Poly-A which furthermore interferes with a photometric quantification of the purified nucleic acids.

## 2.6 Elution procedures

The elution buffer volume does not only influence total yield and concentration of RNA and DNA, but does also influence the ratio between very small and larger oligonucleotides:

### 20 $\mu\text{L}$

The silica membrane is not completely wetted. Only weakly binding very small oligonucleotides like miRNA are eluted efficiently. Larger RNA and DNA are more likely to remain bound to the column. The eluted miRNA is highly concentrated.

### 30 $\mu\text{L}$ (standard)

The standard elution buffer volume of 30  $\mu\text{L}$  is sufficient to wet the silica membrane completely. It results a high total yield of miRNA / RNA / DNA and simultaneously maximizes the concentration.

### 50 $\mu\text{L}$

Increasing the elution buffer volume will further increase the final yield but consequently will reduce the concentration. The gain in yield will usually not compensate for the loss in sensitivity of miRNA detection caused by the dilution of the eluate. Furthermore, larger RNA and DNA will be eluted more efficiently, which might, however, be interesting for the analysis of circulating DNA.



### 3 Storage conditions and preparation of working solutions

#### Attention:

*Buffers MLP and MW1 contain guanidinium thiocyanate. Wear gloves and goggles!*

Storage conditions:

- Store lyophilized RNase-free rDNase at +4 °C on arrival (stable for at least one year).
- All other kit components should be stored at 15–25 °C and are stable until: see package label. Storage at lower temperatures may cause precipitation of salts.

Before starting the first **NucleoSpin® miRNA Plasma** procedure prepare the following:

- **Wash Buffer MW2:** Add the indicated volumes of 96–100 % ethanol to the MW2 concentrate. The buffer can be stored at room temperature (15–25 °C) for at least one year.
- **RNase-free rDNase:** Add the indicated volume of Reaction Buffer for rDNase to the rDNase vial and incubate for 1 min at room temperature. Gently swirl the vial to completely dissolve the rDNase. Be careful not to mix rDNase vigorously as rDNase is sensitive to mechanical agitation. Dispense into aliquots and store at -20 °C. The frozen solution is stable for at least 6 months. Do not freeze/thaw the aliquots more than three times.

NucleoSpin® miRNA Plasma			
REF	10 preps 740981.10	50 preps 740981.50	250 preps 740981.250
Wash Buffer MW2 (Concentrate)	6 mL Add 24 mL ethanol	25 mL Add 100 mL ethanol	100 mL Add 400 mL ethanol
RNase-free rDNase (lyophilized)	1 vial (size C) Add 3 mL Reaction Buffer for rDNase	1 vial (size C) Add 3 mL Reaction Buffer for rDNase	5 vials (size C) Add 3 mL Reaction Buffer for rDNase to each vial

## 4 Safety instructions

When working with the **NucleoSpin® miRNA Plasma** kit wear suitable protective clothing (e.g., lab coat, disposable gloves, and protective goggles). For more information consult the appropriate Material Safety Data Sheets (MSDS available online at [www.mn-net.com/msds](http://www.mn-net.com/msds)).



Caution: Guanidinium thiocyanate in buffer MLP and buffer MW1 can form highly reactive compounds when combined with bleach! Thus, do not add bleach or acidic solutions directly to the sample preparation waste.

The waste generated with the **NucleoSpin® miRNA Plasma** kit has not been tested for residual infectious material. A contamination of the liquid waste with residual infectious material is highly unlikely due to strong denaturing lysis buffer treatment but it cannot be excluded completely. Therefore, liquid waste must be considered infectious and should be handled and discarded according local safety regulations.

### 4.1 Disposal

Dispose hazardous, infectious or biologically contaminated materials in a safe and acceptable manner and in accordance with all local and regulatory requirements.

## 5 Protocol - small RNA and DNA purification from plasma or serum

Before starting with the preparation check that isopropanol is available, that ethanol was added to Wash Buffer MW2, and that rDNase was reconstituted according to section 3.

### 1 Prepare sample


*Note:* See section 2.4 for optional Proteinase K digest.

Add **90  $\mu$ L Buffer MLP** to **300  $\mu$ L sample**.

Vortex for **5 s**.

Incubate for **3 min** at **room temperature**.

*Note:* To process 600  $\mu$ L or 900  $\mu$ L sample material, increase volumes for Buffer MLP, MPP, and isopropanol proportionally. Multiple loading steps will be necessary in step 5. See section 2.3 for more information.

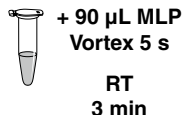


### 2 Precipitate protein

Add **30  $\mu$ L Buffer MPP** and vortex for **5 s**.

Incubate for **1 min** at **room temperature**.

Centrifuge for **3 min** at **11,000 x g** to pellet the protein.

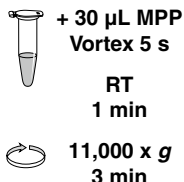

### 3 Transfer supernatant

Transfer the clear supernatant into a new Collection Tube (2 mL, lid).



### 4 Adjust binding conditions

*Note:* Addition of carrier, for example, 2  $\mu$ g of glycogen or 5  $\mu$ g of LPA (linear polyacrylamide), might slightly improve the miRNA yield but usually is not necessary. Poly-A has shown only negligible effects and furthermore might interfere with photometric miRNA quantification.

Add **400  $\mu$ L isopropanol** and vortex for **5 s**.



**5 Bind RNA and DNA**

Place a **NucleoSpin® miRNA Column** in a Collection Tube (2 mL) and load the sample onto the column.

Incubate for **2 min** at **room temperature**.

Centrifuge for **30 s** at **11,000 x g**.

Discard the flowthrough and place the column back into the collection tube.

If more than 300 µL plasma/serum was used, repeat this step until all sample is loaded onto the column.




**6 Optional: DNA digest**

*Note: Co-purified DNA might interfere with qPCR quantification of miRNA. The following on-column digest degrades bound DNA including miRNA genes.*

However, if miRNA specific qPCR detection systems are used or cell free plasma DNA is needed intact for further analysis, skip the rDNase digest and proceed directly with step 7.


**1<sup>st</sup> wash**

Add **700 µL Buffer MW2** to the NucleoSpin® miRNA Column. Centrifuge for **30 s** at **11,000 x g**. Discard flowthrough and place the column back into the collection tube.

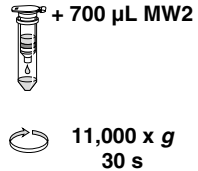

**2<sup>nd</sup> wash**


Add **250 µL Buffer MW2** to the NucleoSpin® miRNA Column. Centrifuge for **2 min** at **11,000 x g**. It is not necessary to discard the flowthrough.

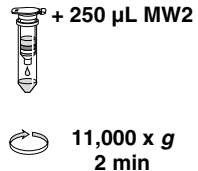

**Digest DNA**

Add **50 µL rDNase** (dissolved in Reaction Buffer for rDNase according to section 3) directly onto the silica membrane of the NucleoSpin® miRNA Column.



Close the lid and incubate at **room temperature** for **15 min**.

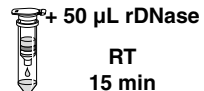

**7 Wash and dry silica membrane****1<sup>st</sup> wash**


Add **100 µL Buffer MW1** to the NucleoSpin® miRNA Column. Centrifuge for **30 s** at **11,000 x g**.

Discard flowthrough and place the column back into the collection tube.

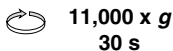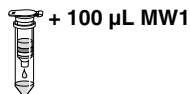**2<sup>nd</sup> wash**

Add **700 µL Buffer MW2** to the NucleoSpin® miRNA Column. Centrifuge for 30 s at **11,000 x g**.

Discard flowthrough and place the column back into the collection tube.


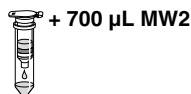
**3<sup>rd</sup> wash**

Add **250 µL Buffer MW2** to the NucleoSpin® miRNA Column.

Centrifuge for **2 min** at **11,000 x g** to dry the membrane completely.

If the liquid in the collection tube has touched the NucleoSpin® miRNA Column after the 3<sup>rd</sup> wash, discard flowthrough and centrifuge again.

*Note: The same collection tube is used throughout the entire washing procedure to reduce plastic waste. If new collection tubes are to be used for each step, see ordering information.*

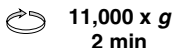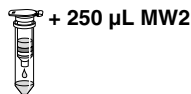**8 Elute RNA and DNA**


Place the NucleoSpin® miRNA Column in a new Collection Tube (1.5 mL).

*Note: The elution buffer volume highly influences the final yield and concentration and, furthermore, influences elution efficiency of large oligonucleotides. See section 2.6 for more information about elution in 20, 30 or 50 µL.*

Add **30 µL RNase-free H<sub>2</sub>O** directly onto the silica membrane of the column.

Incubate for **1 min** at **room temperature**.

Close the lid and centrifuge for **1 min** at **11,000 x g**.

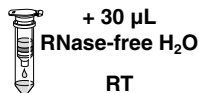

## 6 Appendix

### 6.1 Troubleshooting

Problem	Possible cause and suggestions
Poor or no RNA yield	<i>Reagents not applied or restored properly</i>
	<ul style="list-style-type: none"> <li>• Always dispense exactly the buffer volumes given in the protocols! The correct proportions of buffers MLP, MPP, and isopropanol are essential for optimal yield and purity.</li> <li>• Always follow closely the given instructions with regard to order and mode of mixing (shaking, vortexing etc).</li> <li>• Add the indicated volume of 96 – 100 % ethanol to Buffer MW2 Concentrate and mix thoroughly.</li> <li>• Store kit components at room temperature . Storage at lower temperatures may cause salt precipitation. Heat buffer with precipitated salt to 30 °C until salt is dissolved. Let the buffer cool down to room temperature before use.</li> <li>• Keep bottles tightly closed in order to prevent evaporation or contamination.</li> </ul>
	<i>RNase contamination</i>
	RNA is degraded
Suboptimal performance of RNA in downstream experiments	<i>Inhibition by co-purified RT-PCR inhibitors</i>
	<ul style="list-style-type: none"> <li>• Heme, hemin, and other degradation products of red blood cells strongly inhibit reverse transcription and PCR. Too much plasma or bad plasma quality can result in contamination with these inhibitors. Use less plasma, dilute eluates, perform the optional Proteinase K digest (see section 2.4), or add BSA prior to RT or PCR reactions.</li> </ul>

Problem	Possible cause and suggestions
Suboptimal performance of RNA in downstream experiments (continued)	<i>Carry-over of ethanol or salt</i>
	<ul style="list-style-type: none"> <li>Do not let the flowthrough touch the column outlet after the second MW2 wash. Make sure to centrifuge at the corresponding speed for the respective time in order to remove ethanolic buffer MW2 completely.</li> <li>Check if buffer MW2 has been equilibrated to room temperature before use. Washing at lower temperatures lowers efficiency of salt removal.</li> </ul>
	<p><i>Store isolated RNA properly</i></p> <ul style="list-style-type: none"> <li>Eluted RNA should always be kept on ice for optimal stability since trace contaminations of omnipresent RNases (general lab ware, fingerprints, dust) will degrade the isolated RNA. For short term storage freeze at -20 °C, for long term storage freeze at -70 °C.</li> </ul>

## 6.2 Ordering information

Product	REF	Pack of
NucleoSpin® miRNA Plasma	740981.10/.50/.250	10/50/250 preps
Exosome Precipitation Solution (Serum/Plasma)	740398.12/.20/.60	12/20/60 mL
Exosome Precipitation Solution (Urine)	740399.12/.50/.250	12/50/250 mL
NucleoSpin® miRNA	740971.10/.50/.250	10/50/250 preps
rDNase Set (1 vial rDNase (size F), 7 mL Reaction Buffer for rDNase)	740963	1 set
Proteinase K	740506	100 mg
Collection Tubes (2 mL)	740600	1000

Visit [www.mn-net.com](http://www.mn-net.com) for more detailed product information.

## 6.3 Product use restriction / warranty

**NucleoSpin® miRNA Plasma** kit components are intended, developed, designed, and sold FOR RESEARCH PURPOSES ONLY, except, however, any other function of the product being expressly described in original MACHEREY-NAGEL product leaflets.

MACHEREY-NAGEL products are intended for GENERAL LABORATORY USE ONLY! MACHEREY-NAGEL products are suited for QUALIFIED PERSONNEL ONLY! MACHEREY-NAGEL products shall in any event only be used wearing adequate PROTECTIVE CLOTHING.

For detailed information please refer to the respective Material Safety Data Sheet of the product! MACHEREY-NAGEL products shall exclusively be used in an ADEQUATE TEST ENVIRONMENT. MACHEREY-NAGEL does not assume any responsibility for damages due to improper application of our products in other fields of application. Application on the human body is STRICTLY FORBIDDEN. The respective user is liable for any and all damages resulting from such application.

DNA/RNA/PROTEIN purification products of MACHEREY-NAGEL are suitable for IN-VITRO-USES ONLY!

ONLY MACHEREY-NAGEL products specially labeled as IVD are also suitable for IN-VITRO-diagnostic use. Please pay attention to the package of the product. IN-VITRO-diagnostic products are expressly marked as IVD on the packaging.

IF THERE IS NO IVD SIGN, THE PRODUCT SHALL NOT BE SUITABLE FOR IN-VITRO-DIAGNOSTIC USE!

ALL OTHER PRODUCTS NOT LABELED AS IVD ARE NOT SUITED FOR ANY CLINICAL USE (INCLUDING, BUT NOT LIMITED TO DIAGNOSTIC, THERAPEUTIC AND/OR PROGNOSTIC USE).

No claim or representations is intended for its use to identify any specific organism or for clinical use (included, but not limited to diagnostic, prognostic, therapeutic, or blood banking). It is rather in the responsibility of the user or - in any case of resale of the products - in the responsibility of the reseller to inspect and assure the use of the DNA/RNA/protein purification products of MACHEREY-NAGEL for a well-defined and specific application.

MACHEREY-NAGEL shall only be responsible for the product specifications and the performance range of MN products according to the specifications of in-house quality control, product documentation and marketing material.

This MACHEREY-NAGEL product is shipped with documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. MACHEREY-NAGEL's sole obligation and the customer's sole remedy is limited to replacement of products free of charge in the event products fail to perform as warranted. Supplementary reference is made to the general business terms and conditions of MACHEREY-NAGEL, which are printed on the price list. Please contact us if you wish to get an extra copy.

There is no warranty for and MACHEREY-NAGEL is not liable for damages or defects arising in shipping and handling (transport insurance for customers excluded), or out of accident or improper or abnormal use of this product; defects in products or components not manufactured by MACHEREY-NAGEL, or damages resulting from such non-MACHEREY-NAGEL components or products.

MACHEREY-NAGEL makes no other warranty of any kind whatsoever, and SPECIFICALLY DISCLAIMS AND EXCLUDES ALL OTHER WARRANTIES OF ANY KIND OR NATURE WHATSOEVER, DIRECTLY OR INDIRECTLY, EXPRESS OR IMPLIED, INCLUDING, WITHOUT LIMITATION, AS TO THE SUITABILITY, REPRODUCTIVITY, DURABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, MERCHANTABILITY, CONDITION, OR ANY OTHER MATTER WITH RESPECT TO MACHEREY-NAGEL PRODUCTS.

In no event shall MACHEREY-NAGEL be liable for claims for any other damages, whether direct, indirect, incidental, compensatory, foreseeable, consequential, or special (including but not limited to loss of use, revenue or profit), whether based upon warranty, contract, tort (including negligence) or strict liability arising in connection with the sale or the failure



of MACHEREY-NAGEL products to perform in accordance with the stated specifications. This warranty is exclusive and MACHEREY-NAGEL makes no other warranty expressed or implied.

The warranty provided herein and the data, specifications and descriptions of this MACHEREY-NAGEL product appearing in MACHEREY-NAGEL published catalogues and product literature are MACHEREY-NAGEL's sole representations concerning the product and warranty. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agent or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized; they should not be relied upon by the customer and are not a part of the contract of sale or of this warranty.

Product claims are subject to change. Therefore please contact our Technical Service Team for the most up-to-date information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Applications mentioned in MACHEREY-NAGEL literature are provided for informational purposes only. MACHEREY-NAGEL does not warrant that all applications have been tested in MACHEREY-NAGEL laboratories using MACHEREY-NAGEL products. MACHEREY-NAGEL does not warrant the correctness of any of those applications.

Last updated: 07/2010, Rev. 03

Please contact:

MACHEREY-NAGEL GmbH & Co. KG

Tel.: +49 24 21 969-270

tech-bio@mn-net.com

---


Trademarks / Disclaimer:

NucleoSpin® is a registered trademark of MACHEREY-NAGEL GmbH & Co. KG

TaqMan® is a registered trademark of Roche Molecular Systems, Inc.

All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Regarding these products or services we can not grant any guarantees regarding selection, efficiency, or operation.



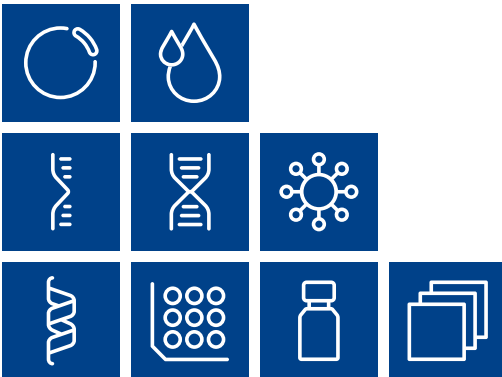

Plasmid DNA

Clean up

RNA

DNA


Viral RNA and DNA

Protein

High throughput

Accessories

Auxiliary tools

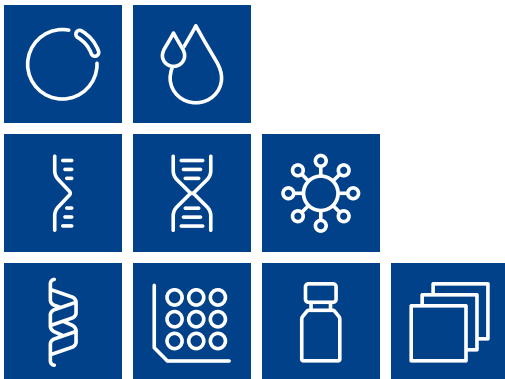

[www.mn-net.com](http://www.mn-net.com)

**MACHEREY-NAGEL**



MACHEREY-NAGEL GmbH & Co. KG  
Valenciennener Str. 11  
52355 Düren · Germany

DE	Tel.: +49 24 21 969-0	info@mn-net.com
CH	Tel.: +41 62 388 55 00	sales-ch@mn-net.com
FR	Tel.: +33 388 68 22 68	sales-fr@mn-net.com
US	Tel.: +1 888 321 62 24	sales-us@mn-net.com



A040554/0420.5