

## Accessories | Syringe repeating adaptor (RAX)

The syringe repeating adaptor (RAX) improves precision and reproducibility when repeatedly dispensing the same volume manually.

## Features and benefits

It has easy to use finger grips and can be set to any desired volume.

The RAX can be used as added protection against both plunger blowout at elevated pressures and plunger bending.

The flat sides allow it to sit securely on the bench with the needle in a raised position helping to prevent damage from rolling.

The RAX is supplied with an innovative double-ended tool to facilitate assembly. It can easily and quickly be installed, and if required it can readily be removed, replaced or swapped between syringes.

## Recommended applications

For repeated manual dispensing of the same volume.


## **Product specifications**

The RAX is suitable for syringes with capacities from 500 nL through to 500  $\mu$ L (except those with plunger stops, or diluter/ dispenser syringes). The one RAX can be used on both 6.5 mm OD and 8 mm OD syringe barrels.

For more information about this product visit www.trajanscimed.com or contact techsupport@trajanscimed.com

Specifications are subject to change without notice.

Analytical syringes manufactured by Trajan Scientific and Medical are intended for analytical and laboratory use only and are not intended or approved for use with food, including the production or packaging of food, nor medical or human in-vivo use.

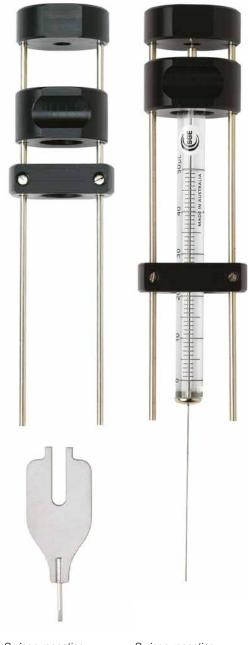

Syringe repeating adaptor (RAX) with assembly tool

Syringe repeating adaptor (RAX) installed on syringe

