


Barcoded well plates

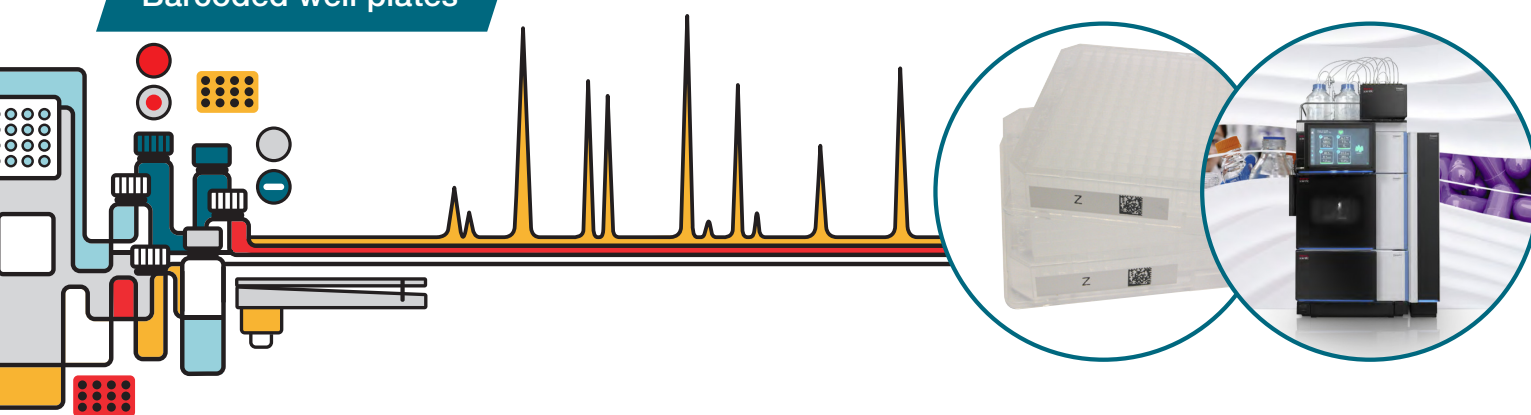

Proactive sample handling solutions

Thermo Scientific™ WebSeal™ well plates, barcoded for Thermo Scientific™ Vanquish™ Horizon and Thermo Scientific™ Vanquish™ Flex UHPLC systems provide advanced sample security and sample integrity. Proactive solutions exist to significantly impact sample security through sample handling solutions such as barcoded well plates for the Vanquish platforms.

For high throughput and/or multiple users with multiple sample formats in an intact sample integrity environment, you need more than fast chromatography. The hardware and consumables need to keep pace with the analysis speed, and our solutions meet that need. The Vanquish Charger can handle up to 20 well plates with an additional three positions in the Vanquish autosampler giving a maximum of 8832 samples. The automation of workflows with barcode reading is also possible.

Thermo Scientific Vanquish Charger features

- The module contains a barcode reader that automatically inventories each plate or rack within the device
- Monitor and manage an inventory of up to nine sample racks, deep well plates, or 20 microtiter plates (15 mm height) in SBS format
- Gain full inventory management via barcode reading with empty shelf level detection and rack/well plate verification

Key features and benefits

Feature	Benefit	Why is it important?
Tested portfolio for demanding HPLC and UHPLC assays	The one and only secure choice for chromatographers	No risk from extractables, no risk of losing valued sample by using an unsuitable plate, which may not be compatible with the sample or solvent
Barcoded well plates for the Thermo Scientific Vanquish platform 2D barcode reader	Assured plate fit with pre-applied orientation barcode in the correct reader position	No risk of placing the plate incorrectly in the autosampler
Secured side orientation in the autosampler	No risk of incorrect side orientation of the plate in the autosampler	No risk of starting the sequence at the wrong plate position, provides sample integrity
Instrument controls the plate orientation and verifies the plate before implementing the sequence	Improved sample security and sample integrity compared to non-barcoded plates	Keeps the sample secure and integrity intact

Why barcoded well plates?

The well established Thermo Scientific WebSeal system offers a new range of chromatography tested, orientation barcoded well plates for the Vanquish Horizon and Vanquish Flex UHPLC systems. The barcodes indicate the correct orientation of the plate to ensure the correct position in the autosampler and charger unit when loading and running the analytical sequence. Plates cannot be accepted in the wrong orientation.

Thermo Scientific WebSeal Well plates barcoded for Vanquish systems (Polypropylene, non-coated, non-sterile)

Description	Total height	Well format	Total volume (µL)	Working volume range µL/well	Pack of	Cat. no.
96-Well Plate, round well, barcoded, polypropylene	31.6	U-shape, 8 mm dia	1300	50-1000	5	60180-P103B
					50	60180-P133B
96-Deep Well Plate, square well, barcoded, polypropylene	44.0	V-shape, square	2000	50-1900	5	60180-P105B
					50	60180-P135B
96-Well Microplate, round well, barcoded, polypropylene	14.7	U-shape, 7 mm dia	270	10-250	100	60180-P206B
					10	60180-P216B
		V-shape, round	220	10-190	100	60180-P207B
					10	60180-P217B
96-Well Plate, round well, barcoded, polypropylene	41.6	U-shape, 7 mm dia	1200	50-900	50	60180-P201B
					5	60180-P211B



Find out more at thermofisher.com/webseal

General Laboratory Equipment—Not For Diagnostic Procedures. ©2021, 2023 Thermo Fisher Scientific, Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. This information is presented as an example of the capabilities of Thermo Fisher Scientific products. It is not intended to encourage use of these products in any manner that might infringe the intellectual property rights of others. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details. **FL21561-EN 0523**